

Graseby 1200 Infusion Pump User Manual

Version: 2.2

MDKMed Medical Technology Co., Ltd.
2025.1.2

Content

1 Symbols, Graphics and Warnings	5
1.1 Descriptions of Graphics and Symbols	5
1.2 Warnings	6
2 Terms and Definitions	12
3 Brief Introduction and Scope of Application	13
3.1 Brief Introduction	13
3.2 Intended Use	13
3.3 Model Naming	14
3.4 Benefits	14
4 Important Features	14
5 Specifications	15
5.1 Basic Specifications	15
5.2 Main Performance	15
5.3 Main Functions and Common Functions	19
6 Structure and Operation Interface	20
6.1 Structural Composition	20
6.2 Display and Operation Interface	22
7 Operation Instructions	22
7.1 Installation of Infusion Pump	23
7.2 Power on and Device Safety Self-test	23
7.2.1 Power on and off	23
7.2.2 Device Safety Self-test	24
7.3 Quick Use Guide	24
7.3.1 Install / Replace IV Infusion Set	24











7.3.2 Select IV Infusion Set	25
7.3.3 Set Infusion Parameters	26
7.3.4 Purge	27
7.3.5 Start Infusion	28
7.3.6 Infusion Completed	29
7.4 Pause or Stop Infusion	30
7.5 Bolus	31
7.5.1 Manual Bolus	31
7.5.2 Hand off Bolus	31
7.6 Lock and Unlock Screen Function	32
7.7 Infusion Mode Selection and Setting	33
7.7.1 Dose Mode	34
7.7.2 Drug Library Mode	35
7.7.3 Drop Speed Mode	36
7.7.4 RTM Mode	37
7.7.5 Sequence Mode	38
7.7.6 Loading Dose Mode	40
7.7.7 Intermittent Mode	41
7.7.8 Micro Mode	41
7.7.9 Feeding Mode	42
7.7.10 Transfusion Mode	44
7.8 View Log	46
8 Alarms	47
8.1 Door Open Alarm	51
8.2 Occlusion Alarm	51
8.3 End of Infusion Alarm	52
8.4 Air-in-line Alarm	52





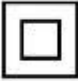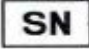

8.5 Battery Empty Alarm	53
8.6 Battery & External Power Disconnect Alarm	54
8.7 Motor Error Alarm	54
8.8 Communication Error Alarm	54
8.9 Battery Error Alarm	55
8.10 KVO Alarm	55
8.11 KVO End Alarm	56
8.12 Standby End Alarm	56
8.13 No AC Power Alarm	56
8.14 Call Back Alarm	57
8.15 Low Battery Alarm	57
8.16 Near End of Infusion Alarm	57
9 System Parameter Setting	58
9.1 Bed Number	58
9.2 Brightness	58
9.3 Alarm Sound Volume	59
9.4 Occlusion Pressure Level	59
9.5 Air Bubble Detection Sensitivity	60
9.6 Bolus Setting	61
9.7 Purge Setting	61
9.8 Call Back Time Setting	61
9.9 KVO Setting	61
9.10 Screen Lock Time	63
9.11 Near End of Infusion Time Setting	63
9.12 Prime Prompt Switch	63
9.13 History Mode Switch	63
9.14 IV Set Brand	64

9.15 Night Mode Setting	64
9.16 Date/Time Setting	65
9.17 Maintenance	66
10 Accuracy Calibration for IV infusion set	66
10.1 Enter the Infusion set Accuracy Calibration Interface	66
10.2 Accuracy Calibration for IV infusion set	66
11 Precautions for Using Disposable Consumables	68
12 Technical Specification	69
13 Restore to Factory Setting	77
14 Use, Maintenance and Removal of the Internal Battery	78
15 Service and Maintenance	79
16 Waste Disposal	80
16.1 Battery	81
16.2 IV infusion set	81
16.3 Infusion Pump	81
17 Electromagnetic Compatibility	81
18 Antistatic Precautions	87
19 Packaging and Accessories	88

1 Symbols, Graphics and Warnings

1.1 Descriptions of Graphics and Symbols

	Caution		Read the User Manual
	Defibrillation prevention Type CF device	RoHS	Compliant to RoHS standards
	Date of manufacturing		Class II device
	Serial Number		Classified collection, uncontrolled discard not allowed
IP24	Ingress Protection Grade		AC (Alternating Current)
	DC (Direct current)		Non-ionizing electromagnetic radiation

	Mute		Lithium battery
	Manufacturer		European Representative
	CE mark demonstrating compliance with RoHS and other EU directives		Medical device marking

1.2 Warnings

Please read the following information carefully, operation that does not strictly follow the guidance will possibly damage the device or do harm to patients' health.

- 1) The infusion pumps are intended to intermittently or continuously delivery the parenteral fluids, medications, blood and blood products via IV infusion route, delivery enteral fluids via the alimentary canal or any route connected to the gastrointestinal system (i.e., the enteral route). The devices can be used together with liquid storage devices /IV infusion sets /Blood transfusion sets /Enteral feeding sets. /Enteral feeding sets.
- 2) Untrained personnel are not allowed to operate the device. The operator must carefully read this User' s Manual, so as to prevent medical accidents caused by improper operation.
- 3) To prevent fire or explosion, it is forbidden to use this device in an environment where flammable or explosive matters are present.
- 4) Do not stack and use other devices that may generate external

radio frequency interference or electromagnetic radiation that may affect the safe operation of this device.

- 5) The operator must use the recommended IV infusion set calibrated in accordance with the requirements described in Section 10 Accuracy Calibration for IV infusion set in this guide, and ensure that the correct IV infusion set brand and type are selected.
- 6) The consumable that is not recommended should never be used for infusion, otherwise it may lead to large infusion inaccuracy and even to become unusable.
- 7) The installation height of this device should not be more than 1 meter above or below the patient's heart.
- 8) It is forbidden to reuse the same IV infusion set on another infusion device.
- 9) This device cannot be used as a portable device.
- 10) It is forbidden to use sharp objects to press on the buttons or the touch screen.
- 11) The Infusion Pump must be serviced and calibrated by trained professional technicians. Before maintenance, make sure to unplug the power cable that supplies power to the device. Untrained personnel are strictly prohibited from opening the device casing, otherwise the eligibility for warranty of the device will be lost.
- 12) Please make sure to use only the parts and accessories provided by MDK.
- 13) When hit hard or dropped, the pump should not be used until it has been checked by trained technical staff.
- 14) Except for wiping the outer surface of the device according to Section 15 Service and Maintenance in this manual, no other part

of the device shall be serviced or maintained by users. If there is any abnormality in the device, please contact the customer service of MDK.

- 15) During the use of this device, the equipment should be placed smoothly and fixedly.
- 16) After loading the IV infusion set, the operator is required to check whether the liquid medicine in the IV infusion set leaks. If there is leakage, stop using the IV infusion set and notify the customer service of MDK.
- 17) Operator should set the infusion parameters strictly based on the doctor's prescriptions. Mistakes in infusion parameter settings may cause harm to patients.
- 18) After setting infusion parameters, the operator must ensure that the infusion device is correctly installed on the device before starting the infusion.
- 19) In order to maintain a high infusion accuracy, the contacting spot of compression on an IV infusion set should be changed every 8 hours.
- 20) To maintain high infusion accuracy, the pump needs to be re-calibrated when there is a significant change in ambient temperature.
- 21) The pump will stop operation automatically when there is an alarm. Press the Start key or click the start button to resume operation after the alarm causing condition is removed.
- 22) To avoid failure or false alarm caused by a dirty occlusion sensor or air-in-line sensor, operator should wipe clean the pump on a regular basis to keep it clean.

- 23) If the sound pressure level of the audible alarm is less than the environmental noise, the operator should turn the alarm volume up to ensure the alarm sound can be heard.
- 24) Pump or accessories may not be usable if their lifetime for use has expired (the lifetime for pump is 8 years). Contact MDK to upgrade to new products.
- 25) The device has a internal rechargeable lithium battery and its lifetime is 2 years.
- 26) Please check the voltage of the internal battery before using it for pump operation. The battery must be replaced and maintained by trained technical personnel in accordance with Section 14 Use, Maintenance and Removal of the Internal Battery in this manual. Replacing the battery by personnel without sufficient training will lead to risks such as over temperature, fire or explosion.
- 27) Please do not connect any other device to the USB and type-c port other than the included DC power adapter shipped with the pump.
- 28) For different types of patients, different occlusion alarm pressure threshold should be set, please refer to the doctor's advice for details.
- 29) Healthcare professional should check on the device during operation on a regular basis, and he/she should also pay attention to medication solution in the IV infusion set before starting the device to make sure the right medicine is in the right infusion channel.
- 30) Please use the roller clamp and other components on the IV infusion set correctly based on the corresponding instruction of the consumable per sec.

- 31) When using this device, please do not plug the power to somewhere that is difficult to plug or unplug. Use an independent power outlet as a measure in case quick disconnection is needed.
- 32) IV infusion set needle is the application part of this product.
- 33) While in normal operation, an alarm will be triggered if the pump door is opened. Please contact MDK for service if this alarm fails to appear.
- 34) If the sticker on the screw hole is removed, then consider the fact that the pump has been tampered with, and discontinue use.
- 35) The product is not AP or APG type device and should not be used in flammable gas environment.
- 36) Don't near active HF SURGICAL device and the RF shielded room of an ME SYSTEM for magnetic resonance imaging, where the intensity of EM DISTURBANCES is high.
- 37) Use of this device adjacent to or stacked with other device should be avoided because it could result in improper operation. If such use is necessary, this device and the other device should be observed to verify that they are operating normally.
- 38) Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this device could result in increased electromagnetic emissions or decreased electromagnetic immunity of this device and result in improper operation.
- 39) Portable RF communications device (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the Graseby 1200, including cables specified by the manufacturer. Otherwise, degradation of the

performance of this device could result.

- 40) The EMISSIONS characteristics of this device make it suitable for use in industrial areas and hospitals (CISPR 11 class A). If it is used in a residential environment (for which CISPR 11 class B is normally required) this device might not offer adequate protection to radio-frequency communication services. The user might need to take mitigation measures, such as relocating or reorienting the device.
- 41) The ME device or ME SYSTEM is suitable for professional healthcare facility environments.
- 42) If the device needs to be used on the move (transport within the hospital): make sure the device is securely fixed and placed. If the device is changed in position, or the pump is severely shaken, the accuracy of the infusion may be affected.
- 43) Do not use unapproved cleaners, materials or chemicals as they may damage device surfaces, labels, or cause device failures.
- 44) Do not route LVP supply bag or administration set right above the pump.
- 45) Do not route the administration set in a way that presents tripping hazard and administration set break off.
- 46) Do not change the height of pump during infusion, otherwise the infusion accuracy may be affected.
- 47) For different types of patients, different occlusion pressure level should be set. For details, please refer to the doctor's advice.
- 48) When the device is powered by the internal battery, the charging indicator light is blue; When the device is powered through the net power supply, the external power indicator light will turn green. At

this time, if the battery is not fully charged, the charging indicator light will turn green at the same time, and the charging indicator light will not turn on when the battery is fully charged.

- 49) Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.
- 50) Do not modify this equipment without authorization of the manufacturer.
- 51) If this equipment is modified, appropriate inspection and testing must be conducted to ensure continued safe use of the equipment.
- 52) If necessary, please contact the company's customer service staff to provide the relevant information of maintenance.

2 Terms and Definitions

Operator: A professionally trained and qualified member of medical staff.

Keep vein open (KVO): After infusion is completed based on the preset parameters, the pump will automatically switch to a mode with extremely low Infusion rate and continue to run (this mode virtually does not have any treatment effect), which is to keep the IV infusion set and vein unobstructed and to avoid the blood flowing backwards.

Intermediate rate: An Infusion rate of 25.00mL/h.

Minimum rate: An Infusion rate of 1.00mL/h.

Free-flow: Drug solution is flowing out in an uncontrolled manner under the effect of gravity.

VTBI: Volume to be infusion.

3 Brief Introduction and Scope of Application

3.1 Brief Introduction

The Graseby 1200 Infusion pump is a high-accuracy infusion device. It is mainly consisted of an electrical control module and a mechanical actuation module, including subsystems such as a control system, a motor driver system, a sensing and monitoring system, an alarm system, a display system, a power system (Adapter) and etc.

3.2 Intended Use

Intended use: The infusion pumps are intended to intermittently or continuously delivery the parenteral fluids, medications, blood and blood products via IV infusion route, delivery enteral fluids via the alimentary canal or any route connected to the gastrointestinal system (i.e., the enteral route). The devices can be used together with liquid storage devices /IV infusion sets /Blood transfusion sets /Enteral feeding sets.


Indication for use: N/A.

Contraindications: None.

Intended patient population: The target population is adults, pediatrics (including newborns) who need intravenous therapy, blood transfusion and enteral nutrition, no other specific requirements.

Intended users of the device: The device is intended to be used by trained healthcare professionals in medical institution environments.

3.3 Model Naming

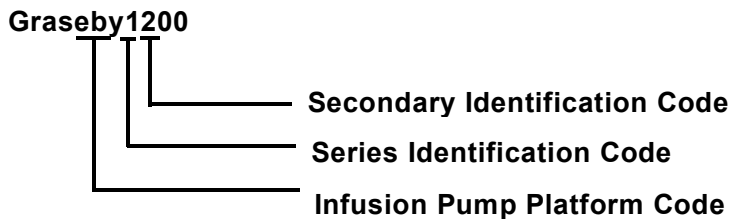

3.4 Benefits

Since the infusion pumps do not come into contact with patient directly, they don't have direct clinical benefit to patients. So the clinical benefit of the infusion pump is that the device enables the precisely administration of intravenous fluid /blood /blood product /enteral fluid procedure to be undertaken, and realizes the monitoring of the infusion processing (the device is able to send alarms during infusion process when occlusion or other circumstances), which can greatly increase the infusion safety.


4 Important Features

- 1) **Accuracy:** The accuracy for infusion rate and volume both are kept within 4.5%.
- 2) **Flow rate:** The Infusion rate can be adjusted from 0.01mL/h to 2000.00mL/h in a continuous manner, which makes the infusion pump capable of meeting various flow rate requirements in different infusion situations.
- 3) **External power supply:** An external power adapter is used, which not only removes the safety concerns of using an internal switching power source but also makes the device lighter and smaller in size.

- 4) **Battery capacity:** The rechargeable internal high-capacity Lithium battery can support normal operation for about 6 hours, which is conveniently helpful during patient transport or power outage.
- 5) **Display:** LCD touch screen display offers high contrast, great visibility and user friendly usability.
- 6) **Occlusion Alarm:** Both upstream and downstream occlusion alarms are available, 10 pressure level is adjustable.
- 7) **Air-in-line alarm:** Based on ultrasonic technology, the device is capable of detecting air bubble sizes down to 25 μ L and initiating air-in-line alarm.

5 Specifications

5.1 Basic Specifications

Dimensions	282 mm×170 mm×120 mm (W×D×H)
Weight	2 kg
Power supply	Network power supply: ~ 100V-240 V, 50/60 Hz Internal battery:  11.1 V rechargeable Lithium battery
Rate of work	40 VA
IV Infusion sets requirements	Refer to Section 11 Precautions for Using Disposable Consumables
Maximum Infusion Rate	2000 mL/h

5.2 Main Performance

Infusion Rate range	0.01 ~2000.00 mL/h with resolution of 0.01 mL/h
---------------------	---

VTBI range	0.01~ 9999.99mL with resolution of 0.01 mL/h
Infusion accuracy	±4.5%
Bolus accuracy	±8.0%
Purge Rate / Bolus Rate	1 mL/h ~ 2000 mL/h, with resolution of 1 mL/h
Purge VTBI / Bolus VTBI	0.10 mL~ 100 mL, with resolution of 0.01 mL/h
KVO Rate	<p>This infusion pump provides two KVO modes to ensure venous patency while guaranteeing infusion safety:</p> <p>1.Constant KVO Mode (V_{KVO} rate setting range: 0.10-5.00 mL/h, minimum step value 0.01 mL/h) -When the infusion rate is $>V_{KVO}$: After the infusion task is completed, the pump automatically runs at the V_{KVO} rate to maintain venous patency. -When the infusion rate is $\leq V_{KVO}$: After the infusion task is completed, the pump continues to run at the infusion rate to maintain venous patency.</p> <p>2.Variable Speed KVO Mode (V_{KVO} rate setting range: 0.10-5.00 mL/h, minimum step value 0.01 mL/h) Users need to set the V_{KVO} rate separately for three scenarios: (a) $V_{KVO}(>10)$ rate for infusion rates > 10 mL/h, (b) $V_{KVO}(1-10)$ rate for infusion rates between 1-10 mL/h, and (c) $V_{KVO}(\leq 1)$ rate for infusion rates ≤ 1 mL/h. The system will automatically matches the KVO</p>

	<p>rate based on the infusion rate:</p> <p>-When the infusion rate is > 10 mL/h: After the infusion task is completed, the pump automatically runs at the $V_{KVO}(>10)$ rate to maintain venous patency.</p> <p>-When the infusion rate is between 1-10 mL/h: After the infusion task is completed, the pump automatically runs at the $V_{KVO}(1-10)$ rate to maintain venous patency.</p> <p>-When the infusion rate is ≤ 1 mL/h and the infusion rate is $> V_{KVO}(\leq 1)$ rate: After the infusion task is completed, the pump continues to run at the $V_{KVO}(\leq 1)$ rate to maintain venous patency.</p> <p>-When the infusion rate is ≤ 1 mL/h and the infusion rate is $\leq V_{KVO}(\leq 1)$ rate: After the infusion task is completed, the pump continues to run at the infusion rate to maintain venous patency.</p>
Infusion Time range	00:00:00 ~ 99:59:59, with resolution of 1 s.
Occlusion threshold	10 levels, with the lowest being $30 \text{ kPa} \pm 20 \text{ kPa}$, and the highest being $120 \text{ kPa} \pm 20 \text{ kPa}$.
Maximum infusion pressure generated by the device	140 kPa
Occlusion alarm	When operated at minimum Infusion rate(1.00

<p>trigger time and Bolus dosage</p>	<p>mL/h): < 1 h when the occlusion alarm pressure threshold is set to the lowest pressure; or < 3h30min when the occlusion alarm pressure threshold is set to the highest pressure.</p> <p>When operated at intermediate rate(25.00 mL/h): < 1min30s when occlusion alarm pressure threshold is set to the lowest pressure, and the Bolus produced during occlusion is < 0.20 mL; < 2min30s, when the occlusion alarm pressure threshold is set to the highest pressure, the Bolus during occlusion is not more than 0.40 mL.</p> <p>(Tested with the Hanaco IV infusion set when an occlusion was created 1 meter away from the pump outlet)</p>
<p>Consumable brand</p>	<p>IV infusion set: 4 brands are recommended, and the default brand is Hanaco, Wego, Yusheng, Shinva. 10 brands can be customized.</p> <p>enteral feeding set: MDK</p> <p>blood transfusion set: Terumo</p>
<p>Supported Infusion modes</p>	<p>11 modes, RVT mode, Dose mode, Drug Library mode, Drop speed mode, RTM mode, Sequence mode, Loading Dose mode, Intermittent mode, Micro mode, Feeding mode, Transfusion mode.</p>
<p>Battery running time</p>	<p>Intermediate rate (25.00mL/h): When fully charged, the battery can run continuously for 6h30min.</p>

	Maximum rate (2000mL/h): When fully charged, the battery can run continuously for 6h10min.
Alarm Mute Time	2 min \pm 10s
Call Back Time	1 min ~ 60 min \pm 10s
Classification	Type II CF continuous operating volumetric Infusion pump with internal power supply; Grade IP24, non AP/APG type device.
Ambient temperature and humidity	Ambient temperature of transportation and storage: -20 °C ~ + 55 °C Ambient temperature for operation: +5 °C ~ + 40 °C Ambient humidity for transportation, storage and operation: 20% ~ 90% Ambient pressure for transportation, storage and operation: 700 hPa ~ 1060 hPa
Software version	Graseby 1200-V1
Service lifetime	8 years

5.3 Main Functions and Common Functions

- 1) Set infusion rate, infusion VTBI and real-time data display function;
- 2) Display of completed infusion volume;
- 3) Purge / Bolus;
- 4) Alarms;
- 5) The Infusion rate will be automatically changed to KVO Rate after the VTBI complete alarm is activated
- 6) Temporary mute for alarm sound and timer for alarm sound recovery;
- 7) Automatic free-flow stopping function;

- 8) Displays the accumulated quantity infusion and supports clearance
- 9) A variety of brands for IV infusion set are supported;
- 10) Built-in battery;
- 11) External power adapter;
- 12) Closed-loop infusion protection safety function
- 13) Contains Dose-Error Reduction Software

6 Structure and Operation Interface

6.1 Structural Composition

The structure of the device consists of a pump casing, a motor drive system, an input system, a storage system, a control system, a display system, a sensing and monitoring system, an alarm system and a power supply system.

Names for parts and components:

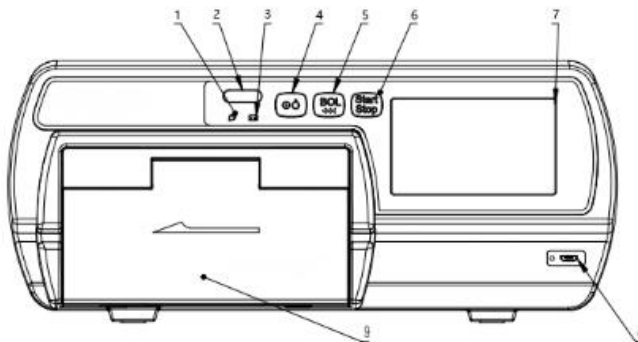


Figure 6-1-1 Front view

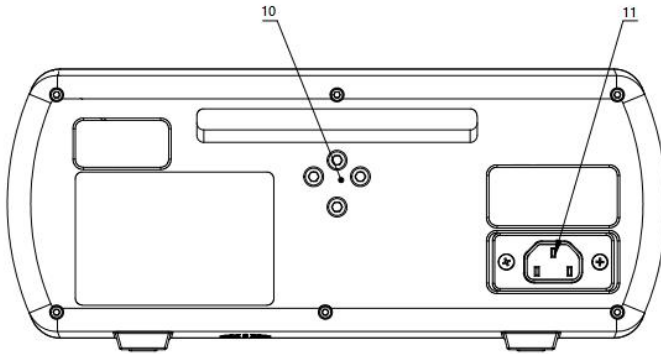


Figure 6-1-2 Back view

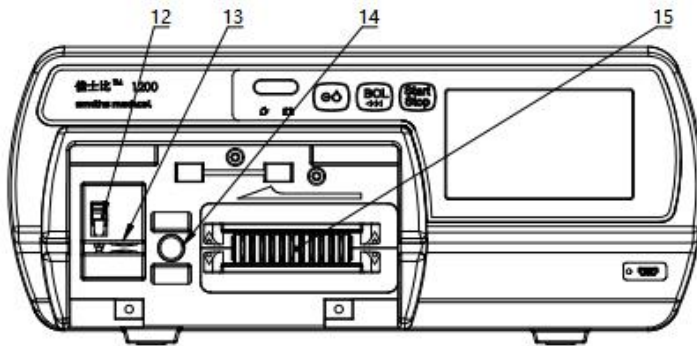


Figure 6-1-3 Pump door open diagram

1	External power indicator	2	Operation status indicator	3	Charging indicator
4	Power on/off key	5	Purge/bolus key	6	Start/stop key
7	Touch screen display	8	Drop sensor port	9	Pump door lock handle
10	Fixed clamp screw hole	11	Power Port	12	Liquid stop clamp

13	Air in line sensor	14	Occlusion sensor	15	Pump disc
----	--------------------	----	------------------	----	-----------

6.2 Display and Operation Interface

The display interface is shown below.

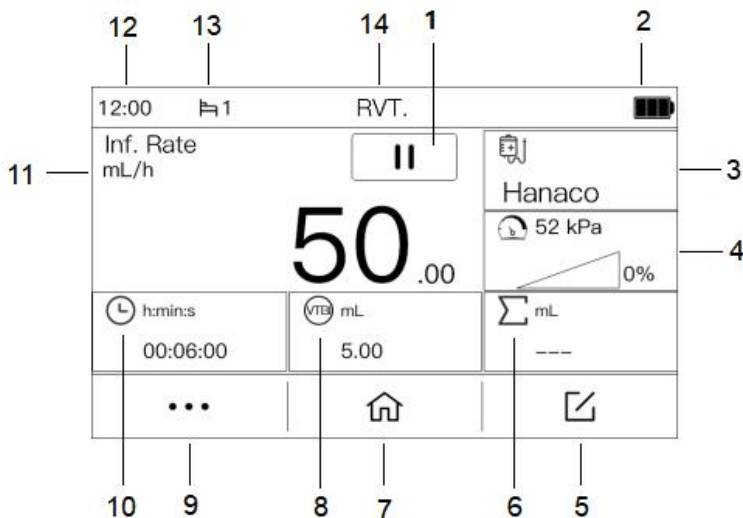


Figure 6-2 Operation interface on the screen

1	Infusion status	2	Battery	3	Brand
4	Occlusion pressure level	5	Edit	6	Total infused volume
7	Home	8	Volume to be infused	9	More
10	Time remaining	11	Infusion rate	12	Time
13	Bed No.	14	Infusion mode		

7 Operation Instructions

Install Infusion pump → Power on → Device safety self-test
 → Install IV infusion set → Select IV infusion set brand →
 Parameters setting → Prime / Purge → Start infusion → Infusion
 completed → Remove accessories → Power off.

Before infusion starts, please confirm that the IV infusion set in use matches the current IV infusion set setting selected in the menu. Any IV infusion set which brand is not included in the list of recommended brands must be calibrated before being used.

7.1 Installation of Infusion Pump

First loosen the locking screw, install the clamp to the pole of the infusion stand, adjust the height of the clamp, and then tighten the locking screw. The operator must make sure that the Infusion pump is positioned in a secure, stable and reliable manner.

7.2 Power on and Device Safety Self-test

7.2.1 Power on and off

Before connecting to the mains power, check if there is any foreign matter inside the power outlets (such as drug solution residue).

Connect to the mains power, check the power indicator on the pump front panel. If the indicator is not lit up, check the connection of power cable and the pump, or check if there is a power outage. Then press the power Key on the front panel to turn the device on.

After infusion therapy is completed, press the power key and click Power-Off button to turn the device off. Do not power off when the device is in operation mode, otherwise the infusion therapy will

be stopped.

7.2.2 Device Safety Self-test

The pump will perform an automatic safety self-test after powered on, if the test is passed then there will be two short beeps and the operation status indicator will be lit up in stable green color. If a continuous alarming sound is initiated or there is no any sound at all, then the device cannot be used, please contact the customer service immediately.

7.3 Quick Use Guide

7.3.1 Install / Replace IV Infusion Set

Firstly, hold the pump door lock handle from the opening below the pump door and pull it outward. Open the pump door and push the liquid stop clamp upwards to make it open. After straightening the infusion tube under the drip pot of the infusion set, horizontally clip it into the gap at both ends of the infusion pump, straighten the infusion set, and place it in the positioning groove (infusion direction should be indicated according to the infusion direction). Lift the pump door lock handle upwards to ensure that the hook has fastened the door pin, Press down on the lock handle again to close the pump door, and the surface of the lock handle should be flush with the pump casing. Rotate the liquid stop pulley on the infusion set to the open state, and the installation of the infusion set is completed.

Before changing the clamping position or medication of the infusion device, the liquid stop pulley on the infusion device must be

turned to the closed state to prevent the occurrence of drug self flow. To replace or reinstall an infusion device, first open the pump door, push up the liquid stop clamp to keep it open, remove the infusion device, reinstall the infusion device, and then rotate the liquid stop pulley on the infusion device to open it.

7.3.2 Select IV Infusion Set

After the infusion pump is powered on and the safety self-test is passed, the parameter setting page will show up. Click the Brand button in the upper right corner to enter the IV infusion set brand selection page.

After clicking an IV infusion set brand to make a selection, the system will automatically return to the parameter setting page.

Please check if the IV infusion set displayed on the right side of the screen matches the set that is being used.

It is possible that the IV infusion set from the same brand may have different characteristics if they are from different lots, which will affect their infusion accuracy if they are not calibrated before use. In that case, calibration of the IV infusion set is recommended, which is described in Section 10.2 Accuracy Calibration for IV infusion set.

IV Set Brands	
Hanaco	
< > ↶	

Figure 7-3-2 Select IV infusion set

7.3.3 Set Infusion Parameters

General method:

When the infusion pump is standby, click "☑" on the touch screen to enter the RVT mode parameter setting interface. Click "Inf. rate" on the touch screen, a numeric button board appears on the screen, click to enter the value of the Infusion rate to be set, and press "√" on the screen to complete the input.

Setting the VTBI and infusion time is the same as setting the infusion rate above. After all parameters are set, click the "√" button to confirm the parameters.

Quick setting method:

When the infusion pump is standby, click "Inf. rate" value on the screen, and a numeric button board appears. Click to enter the value of the Infusion rate to be set, and press "√" on the screen to complete the input.

RVT.		
Inf. Rate	---	mL/h
VTBI	---	mL
Inf. Time	--:--:--	H:M:S
✓	<	> ↶

Figure 7-3-3-1 Set infusion parameters

1	2	3	⊗
4	5	6	↶
7	8	9	✓
.	c	0	

Figure 7-3-3-2 Keyboard input value

7.3.4 Purge

Manual purge:

After the infusion parameters are set, confirm that the infusion line is disconnected from the patient and long press the Bolus button. The device will operate at the default purge rate to quickly empty the air in the infusion line. Release the Bolus button to exit the purge state.

Auto purge:

After the infusion parameters are set, confirm that the infusion line is disconnected from the patient and press the Bolus button.

The device will pop up a prompt to confirm that the infusion line is disconnected and request for purge. After press " ✓ " on the screen, the device will run at the default purge rate and default total volume set by the system to quickly empty the air in the infusion line. Press the stop button again and the purge will be stopped.

The purge volume is not included in the Infusion accumulation. When the Purge is running, the Air-in-line alarm is not suppressed, and the other alarms are normal.

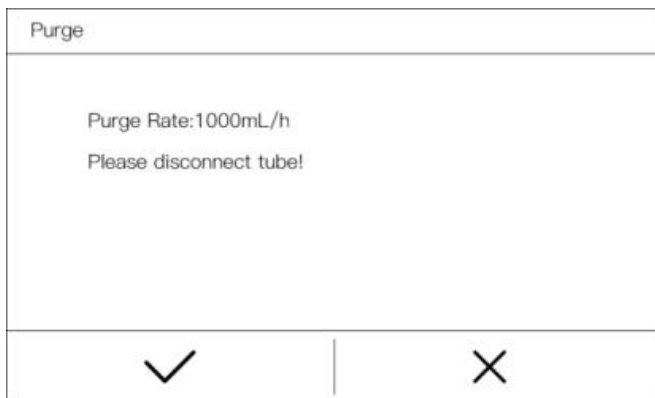


Figure 7-3-4 Confirm to purge

7.3.5 Start Infusion

Press the Start/Stop key or click the start button and the pump will start to run according to the set infusion parameters, as shown in the following figure.

When the device is working, after revising the infusion flow rate in real time, the motor will synchronously change the voltage to increase the motor speed, so that the pump can synchronously reach the changed flow rate.

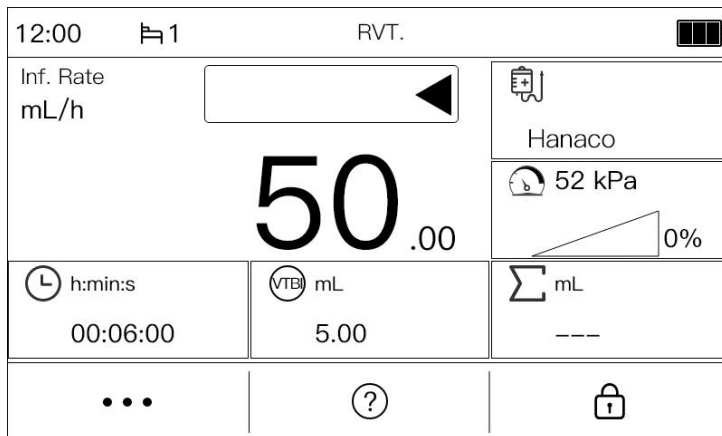


Figure 7-3-5 Infusion interface

7.3.6 Infusion Completed

The infusion is completed when the infusion accumulation volume reaches the VTBI set for the infusion task.

If KVO is enabled, the device will convert to the KVO Rate to continue running automatically and trigger the "KVO" high priority alarm at the same time. Click "√" on the screen to exit the KVO infusion state.

If KVO is not enabled, the device triggers the "Infusion complete" high priority alarm with an alarm sound. Click the "√" on the screen to remove the high priority alarm.

After the infusion task is completed, remove the infusion attachment, follow the method of replacing the infusion set, remove

the infusion set that is no longer in use, press the "Power" button, and click "Shutdown" to turn off the device. Pull the circular ring on the fixed clamp base upwards and remove the infusion pump by pulling it outward.

7.4 Pause or Stop Infusion

Normal infusion state is as shown in Figure 7-3-5.

Press the Start/Stop key or click the stop button during infusion operation can pause the operation, as shown in the following figure.

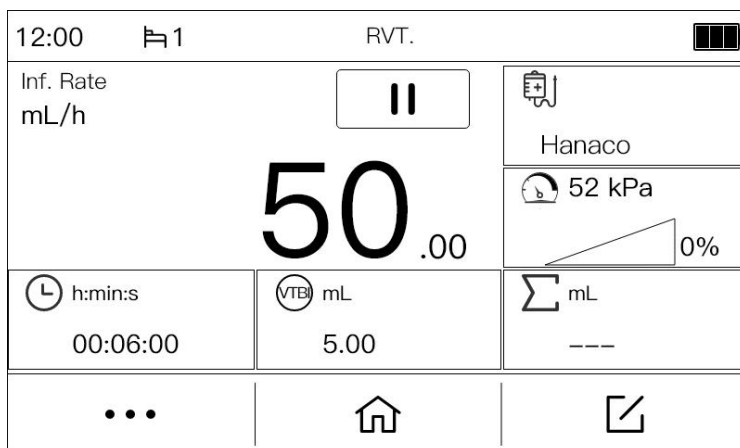


Figure 7-4 Pause or Stop Infusion

On the Infusion pause page, press the Start/Stop key or click the start button again can start the device operation, and the device will continue to infuse the remaining volume based on the set Infusion rate.

During the infusion pause process, if modify any parameters such as infusion rate, total amount, and time is considered as a new infusion task. When the start stop button is pressed again, the

device will start the new infusion task according to the new infusion parameters.

When the device triggers an alarm, it emits an alarm sound. Click the "Alarm Pause" button on the screen to pause the alarm. After 2 minutes, if the alarm source is not released, the alarm sound will automatically resume.

When want the device to stop infusion, press the start stop button to stop.

7.5 Bolus

7.5.1 Manual Bolus

In the infusion state, long press the Bolus key, the infusion pump enters bolus state, and runs according to the rate set in "Setting - Bolus - Bolus Rate". Release the bolus key to stop the bolus infusion, and the infusion pump returns back to the normal infusion state before the bolus to continue infusion, the bolus volume is included in the infusion accumulative volume.

7.5.2 Hand off Bolus

In the infusion state, press the Bolus key or click the "BOL" button, enter the Hand Off Bolus page, set the bolus parameters, click the "✓" button, the infusion pump enter into bolus state until the bolus VTBI is completed, the infusion pump returns to the normal infusion state continue the infusion, the bolus volume is included in the infusion accumulative volume.

Hand Off Bolus		
Bolus Rate	---	mL/h
Bolus VTBI	---	mL
Bolus Time	---:---:---	H:M:S
✓		↶

Figure 7-5-2-1 Bolus setting interface

12:00	1	RVT.	
Bolus Rate mL/h	<input type="text" value="2000.00"/>		
h:min:s 00:00:09	mL 5.00	mL ---	
...	?		

Figure 7-5-2-2 Bolus infusion interface

7.6 Lock and Unlock Screen Function

The device automatically locks the screen after running for a period of time. When the device is in the lock screen, click the screen and a prompt will pop up asking whether to unlock the screen, click the " ✓ " button to unlock the screen.

Auto Lock time settings See Home - Setting - Auto Screen Lock.

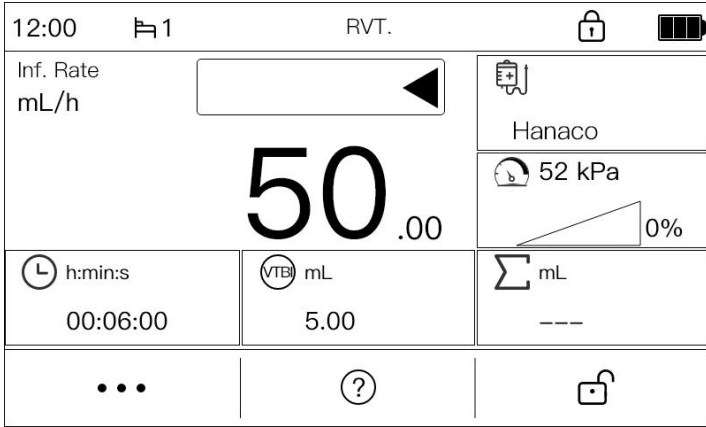


Figure 7-6 Lock screen

7.7 Infusion Mode Selection and Setting

Inf. Mode	1/2
RVT.	Drug Library
Loading Dose	Micro
Dose	RTM
Sequence	Intermittent
< > ↶	

Figure 7-7 Infusion Mode Selection and Setting

Except for the RVT mode on the home screen, there are 10 infusion modes on the Infusion Mode page: Dose mode, Drug Library mode, Drop speed mode, RTM mode, Sequence mode, Loading Dose mode, Intermittent mode, Micro mode, Feeding mode, Transfusion mode.

In the RVT setting page, the infusion rate, VTBI and infusion time can be set in a variety of combinations, forming the following four combinations of infusion mode: rate + volume (R+V) Mode, rate + time (R+T) Mode, volume + time (V+T) Mode, rate (R) Mode. Therefore, there are 14 different infusion modes for the device in total.

The setting of the RVT mode should follow the instructions in the 7.3 "Quick Use Guide" above. The settings for the other modes are outlined below.

7.7.1 Dose Mode

Enter the Dose mode settings interface, as shown in Figure 7-7-1. After setting the Dose, Solution, Concentration, Dose Rate, Weight, VTBI and Infusion Rate, click " ✓ " to confirm the parameters.

In the main interface of Dose mode, click the "unit" in the upper left corner to switch the display of "Dose Rate" and "Infusion rate".

Dose		1/2
Dose	---	mg
Solution	---	mL
Conc.	---	mg/mL
Dose Rate	---	mg/kg/h

Dose		2/2
Weight	---	kg
VTBI	---	mL
Inf.Rate	---	mL/h

✓		<		>		↶
---	--	---	--	---	--	---

Figure 7-7-1 Dose mode setting

In the Dose Mode Settings interface, click "Unit" to the right of any parameter in "Dose", "Concentration" and "Dose Rate" to select different unit expression modes, and the other two corresponding units will be automatically adjusted.

Concentration and Infusion Rate are calculated as follows:

Concentration calculation formula:

$$\text{Concentration(mg/mL)} = \frac{\text{Dose (mg)}}{\text{Solution (mL)}}$$

Infusion rate calculation formula:

Infusion Rate(mL/h)


$$= \frac{\text{Dose Rate (mg/kg/h)} \times \text{weight (kg)} \times \text{Solution (ml)}}{\text{Dose (mg)}}$$

$$\text{Infusion Rate(mL/h)} = \frac{\text{Dose Rate (mg/h)} \times \text{Solution (ml)}}{\text{Dose (mg)}}$$

7.7.2 Drug Library Mode

Enter the drug library mode settings interface. Figure 7-7-2.

Select the name and the specific specifications of the drug that requires infusion. The device will enter the drug library mode settings page and automatically brings in the drug-related parameters. At this point, the drug name and drug specifications are displayed in the title bar of the parameter settings page. After setting the parameters, click " ✓ " to confirm the parameters.

ADREnaline		1/2
Conc.	---	mg/mL
Weight	---	Kg
Dose Rate	---	mg/Kg/mL
Inf. Rate	---	mL/h
		


ADREnaline		2/2
VTBI	---	mL
		

Figure 7-7-2 Drug library mode setting

The Drug library mode has built-in DERS (Dose-error Reduction Software) functionality to reduce medication errors and improve infusion safety.

7.7.3 Drop Speed Mode

Enter the system maintenance page, and enable the drip transfer function.

The drop speed mode setting page is shown in Figure 7-7-3. After the drop speed parameter is set, the system will automatically convert it to a flow rate and display it. The other steps for the settings in this mode are the same as those in the RVT mode described in the Quick Start Guide section.

Drip		
VTBI	---	mL
Drop Rate	---	dot/min
Inf. Rate	---	mL/h
✓		↶

Figure 7-7-3 Drop speed mode setting

7.7.4 RTM Mode

Enter the RTM mode setting interface, as shown in Figure 7-7-4.

RTM.	1/2					
Total Time	--:--:--	H:M:S				
Up Time	--:--:--	H:M:S				
Down Time	---:--	H:M:S				
VTBI	---	mL				
✓		<		>		↶

RTM.	2/2		
Plateau Rate	--- mL/h		
✓	<	>	↶

Figure 7-7-4 RTM mode setting

In the RTM mode parameters, VTBI, Up Time, and Down Time must be set. After setting one of the two parameters of Plateau rate or Total Time, the other parameter will be calculated automatically. After setting the parameters, click "✓" to confirm the parameters.

When the device is started, the infusion rate gradually increases from 0 to Plateau Rate during the Up Time and then maintains the rate. When the remaining time is equal to the Down Time, the infusion rate gradually decreases until it reaches 0 and the infusion is completed.

7.7.5 Sequence Mode

Enter the Sequence Mode settings interface to set the number of sequences (up to 10 groups), click the ">" button to enter the infusion parameter settings interface, as shown below.

Sequence		1/3	
Seq. Count	1		
✓	<	>	↶

Sequence		2/3	
S1-Inf. VTBI	---	mL	
S1-Inf. Time	--:--:--	H:M:S	
S1-Inf. Rate	---	mL/h	
✓	<	>	↶

Sequence		3/3	
Total VTBI	---	mL	
Inf. Time	--:--:--	H:M:S	
✓	<	>	↶

Figure 7-7-5 Sequence mode setting


As shown in the figure, Set any two of the Infusion VTBI, Infusion Rate, Infusion Time, the device will automatically calculate another parameter.

Click the ">" button at the right bottom, set 1 to 9 different infusion parameters according to clinical needs, and click the "✓" button to confirm the infusion parameters.

hen the device completes the infusion parameters of the first sequence, it automatically switched to the parameter run of the second sequence until the set parameters of all sequences were completed and the infusion is completed.

7.7.6 Loading Dose Mode

Enter the Loading Dose mode settings interface. as shown in Figure 7-7-6.

Loading Dose		1/2
VTBI	---	mL
Loading VTBI	---	mL
Loading Rate	---	mL/h
Maintain Rate	---	mL/h
		


Loading Dose		2/2
Loading Time	--:--:--	H:M:S
Maintain Time	--:--:--	H:M:S
		

Figure 7-7-6 Loading Dose Mode setting

As shown in the figure, After setting the VTBI, Loading VTBI, Loading Rate, and Maintain Rate, the device automatically calculated the Loading Time and Maintain Time, and click the " ✓ " button to confirm the infusion parameters.

7.7.7 Intermittent Mode

Enter intermittent mode settings interface.

As shown in the figure, set the Single VTBI, Single Rate, Intermittent Time and Maintain Rate, and click the " ✓ " button to confirm the infusion parameters. After the device is started, the infusion will start at the Single Rate. When the infused volume is equal to the single VTBI, the device will automatically continue infusion according to the Maintain Rate. When the Maintain Rate is set to 0, the device will run at 0mL/h. The device runs at Maintain Rate until the time is equal to the Intermittent Time, the device automatically switches to the Single Rate to continue the infusion, so as to cycle.

Intermittent		
Single VTBI	---	mL
Single Rate	---	mL
Inter Time	--:--:--	H:M:S
Maintain Rate	---	mL/h
✓		↶

Figure 7-7-7 Intermittent mode setting

7.7.8 Micro Mode

Enter Micro mode settings interface.

Set any two of the Infusion Rate, VTBI, and Infusion Time, the device will automatically calculate another parameter. The Infusion Rate should not exceed 100mL/h. After the device is started, the infusion starts at the Infusion Rate. When the infused volume is

equal to the VTBI, the infusion will be stopped automatically.

When only sets the Infusion Rate, the device runs at the Infusion Rate until the operator stops the infusion or the device triggers a high priority alarm to stop the infusion.

7.7.9 Feeding Mode

Enteral nutrition infusion refers to the technique of precisely controlling the administration of enteral nutrition fluids (such as intact protein formulas, short peptide formulas, etc.) into a patient's gastrointestinal tract via a nasogastric tube, nasoenteric tube, or stoma tube using an infusion pump. This method is suitable for patients with dysphagia, normal gastrointestinal function but unable to eat independently, to maintain their nutritional and metabolic needs. Enteral nutrition fluids must be infused at a constant rate to avoid gastrointestinal intolerance (such as diarrhea, bloating) or the risk of aspiration. The infusion pump can accurately regulate the rate, especially for critically ill patients, postoperative patients, or those requiring long-term nutritional support, ensuring the efficiency and safety of nutrient absorption.

Mode Features:

Precise Rate Control: Supports fine adjustments down to 0.1 mL/h to meet varying flow rate requirements from the initial adaptation period to the stable phase.

Safety Alerts: Equipped with alarms for tube blockage, completion of infusion, and abnormal pressure to reduce the risk of tube dislodgement.

Wide Compatibility: Compatible with various viscosities of nutrition fluids and standard enteral feeding tubes, reducing issues of sedimentation or blockage.

Precautions:

- Confirm the tube position (e.g., in the stomach or intestine) before infusion to avoid complications from misplacement.

- Shake the nutrition fluid well and control its temperature (recommended to be close to body temperature). Rinse the tube with warm water before and after infusion to prevent blockage.

- Start the initial infusion at a low rate (e.g., 20-50 mL/h) and gradually increase it, monitoring the patient's tolerance.

- For long-term use, regularly replace the tube and strictly follow aseptic procedures to avoid contamination.

Install the special nutrition tube according to the operation instructions in Section 7.3.1, and select the consumables as the special nutrition tube according to the operation instructions in Section 7.3.2. Click Home-Infusion Mode-Feeding Mode to enter the page for mode parameter settings. The total infusion volume (mL), flow rate (mL/h), single infusion duration, and pressure alarm threshold need to be set. It is recommended to adjust the flow rate in stages until the target rate is reached.

According to the method in Section 7.3.3, set two out of the three items: feeding flow rate, feeding quantity, and feeding time, and the device will automatically calculate the third item. When the feeding flow rate is set and both the feeding quantity and feeding time are set to 0, the set speed will be maintained until the fluid bag

is empty. After setting the feeding parameters and removing air bubbles from the infusion line according to Section 7.3.4, press the Start-Stop button or click the Start button to begin the infusion.

7.7.10 Transfusion Mode

The transfusion mode refers to the process of precisely controlling the infusion rate and volume of blood or blood components (such as red blood cells, plasma, platelets, etc.) using an infusion pump to ensure the safe and efficient administration of blood products to patients.

Blood transfusion requires strict control of the flow rate and total volume. Infusing too quickly may lead to risks such as circulatory overload and hemolytic reactions, while infusing too slowly may affect the timeliness of treatment. The infusion pump achieves constant-rate infusion through preset parameters, especially for critically ill patients, children, and scenarios requiring precise dose control.

Mode Features:

Precise Rate Control: Supports fine adjustments down to 0.1 mL/h to meet the needs of different patients.

Intelligent Protection: Equipped with air bubble detection, tube blockage alarms, and pressure limit alarms to prevent transfusion accidents.

Convenient Operation: Features a built-in blood transfusion-specific mode to simplify parameter settings and improve emergency efficiency.

Precautions:

- Strictly verify blood type, expiration date, and patient information before transfusion.

- Use the recommended blood transfusion-specific tubing to avoid hemolysis caused by material issues with regular infusion tubing.

- Blood products should be warmed to room temperature in advance, and vigorous shaking is prohibited.

- Monitor the patient's vital signs throughout the process and be vigilant for allergic reactions or transfusion-related complications.

Install the special blood transfusion tube according to the operation instructions in Section 7.3.1, and select the consumables as the special blood transfusion tube according to the operation instructions in Section 7.3.2. Click Home-Infusion Mode-Transfusion Mode to enter the mode parameter setting page. During operation, the following must be clarified: total infusion volume (mL), flow rate (mL/h), and pressure alarm threshold (adjusted according to vascular conditions). It is recommended to start with a low flow rate and gradually adjust to the target rate after observing no adverse reactions.

The method for setting transfusion mode parameters is the same as that for the RVT mode. The three parameters of flow rate, infusion volume, and infusion time can be set to form the following four combination modes: flow rate + total volume mode, flow rate + time mode, total volume + time mode, and RVT mode (without total volume and time).

After setting the transfusion parameters and removing air bubbles from the infusion line, press the Start-Stop button or click the Start button to begin the infusion.

7.8 View Log

On the Home - Event Log Page, event logs such as device infusion status and alarm can be displayed. Click this event can view the detailed event information such as Infusion Rate, VTBI, time, Alarm priority and time.

When the pump log store reaches the upper limit of the pump capacity, the oldest log will be overwritten by the new log.

When the alarm system is powered off, the log still exists.


EventLog	1/2
Low Battery	2020-01-01 >
Call Back Alarm	2020-01-01 >
Low Battery	2020-01-01 >
Low Battery	2020-01-01 >
	

Figure 7-8-1 Event log


No AC Power	
Time	10-05 21:46:16
Alarm priority	Low priority alarm
	

Figure 7-8-2 Alarm priority

8 Alarms

Alarm refers to the infusion changes caused by the abnormal infusion circuit or the failure of the Infusion pump itself, which leads to the failure of the infusion to the patient. The Infusion pump prompts the medical staff through sound, light, screen signs and other ways.

Alarm classification prompts of the device:

No.	Alarm	Priority	Alarm category	Alarm conditions
8.1	Door Open	High	Latching	The pump door is not closed during operation or purge.
8.2	OCCL (Occlusion)	High	Latching	When the infusion line is occluded.
8.3	End of Infusion	High	Latching	When the infused volume reaches to the VTBI.
8.4	Air-in-line	High	Latching	Air bubbles are detected in the line.
8.5	Battery Empty	High	Latching	When the internal battery is running out.
8.6	Battery&External Power Disconnect	High	Unlatching	When the device is running, the battery and external power is disconnected at the same time.
8.7	Motor Err.	High	Latching	In the event of a motor failure.

8.8	Com. Err. (Communication error)	High	Latching	Monitor the CPU for communication handshake errors.
8.9	Battery Error	High	Latching	The device does not detect battery signal or battery disconnect when plug in the external power.
8.10	KVO	High	Latching	The infusion is complete with KVO is enabled.
8.11	KVO End	High	Latching	KVO status run for 30 minutes until the KVO task is complete.
8.12	Standby End	High	Latching	Standby End
8.13	No AC Power	Low	Unlatching	When the device is disconnected from the external power and operated with batteries.
8.14	Call Back	Low	Unlatching	When the device is on and not running, no operation is performed within the set time
8.15	Low Battery	Low	Unlatching	When the internal battery power is lower than 30%.

8.16	Near End of Infusion	Low	Unlatching	When the remaining time is less than or equal to the set near end of infusion time.
------	----------------------	-----	------------	---

The device alarm indicator characteristics:

Alarm priority	Indicator color	Flicker frequency		Rate
High priority	Red	2 Hz	0.7 Hz (Battery & External Power Disconnect)	50%
Medium priority	/	/		/
Low priority	Yellow	Normally turned on		100%

High-priority and low-priority alarms are distinguished in sound and light according to standard requirements. When an alarm occurs, the operator can accurately detect it at 1 m away from the alarm system. The delay time of triggering the alarm signal is not more than 2 s.

After powering on, the status indicator lights up, and the device automatically conducts a safety check. After passing the self-check, you will hear two short beeps of "DiDi", which means that the alarm system is normal. If the status indicator does not light up or you hear a continuous alarm sound or no prompt sound after booting, it means that the alarm system is faulty and the device cannot be used normally. It can be put into use after being repaired.

The sound pressure range of the audible alarm signal is 60-95dB.

Note that this device prohibits access to the change or storage

change alarm function. In the process of adjusting the alarm limit or alarm preset, the operation of the alarm system still runs according to the last setting. This device alarm is a technical alarm state.

=====



Caution

Latching alarm signal: The alarm signal that continues to be generated after the trigger event no longer exists, and does not stop until the operator deliberately acts (click the "√" button);

Unlatching alarm signal: When the related trigger event is no longer When it exists, automatically stop the alarm signal generated.

device alarm announcement sequence:

High priority alarm sound priority principle, that is, when the device is in the low priority alarm sound state, when a high priority alarm is generated, the original low priority alarm sound is interrupted, the high priority alarm sound is broadcast, and the high priority alarm sound is displayed at the same time Level alarm prompt information.

The device is in a high priority alarm. When a low priority alarm is generated, the high priority alarm continues to broadcast without being interrupted.

The device is in low priority alarm. When a low priority alarm is generated, it still reports a low priority tone and displays the latest alarm prompt information.

The device detects that a visual alarm and an audible alarm appear immediately.

When the power loss time is less than 30 seconds, the alarm settings before the power loss will automatically restore.

=====

8.1 Door Open Alarm

Cause: When the infusion pump is running, if the pump door is not closed, or the pump door is opened by accident, the device will trigger an alarm, stop running, make a high priority alarm sound, the screen appear with the message "Door Open" and the operation status indicator flash red at the same time.

Solution: Click the "✓" on the screen to clear the alarm, the word " Door Open " disappears, and returns to the infusion pause interface. Check the pump, close pump door and continue to operate.



Figure 8-1 Door open alarm

8.2 Occlusion Alarm

Cause: When the infusion line is occluded, occlusion sensor detects that it is exceeding the set value, the device will trigger an alarm, stop running, make a high priority alarm sound, the screen appear with the message "OCCL" and the operation status indicator flash red at the same time.

At the same time, as a infusion safety protection mechanism, the motor reverses back to pump a small amount of liquid medicine

to reduce the dose of the bolus before occlusion relief.

Solution:

- 1) Click the " ✓ " button on the screen to clear the alarm and the message " OCCL " disappear.
- 2) Check whether the IV infusion set line is kinked, whether the patient presses into the infusion line and other issues, eliminate the problem and restart the infusion.
- 3) If there is still an occlusion alarm, shut off the roller clamp on the IV infusion set, open the pump door, pull out the IV infusion set, check whether the filter or the needle on IV infusion set is occluded, change to a new IV infusion set if necessary and restart infusion.

8.3 End of Infusion Alarm

Cause: If KVO is disabled, When the infused volume reaches the VTBI, the device will trigger an alarm, stop running, make a high priority alarm sound, the screen appear with the message "End Of Infusion" and the operation status indicator flash red at the same time.

As a infusion security mechanism, If KVO is enabled, the device will automatically convert to KVO Rate to continue the infusion.

Solution: Click the " ✓ " button on the screen to clear the alarm and the message "End of Infusion" disappears. The device can be set up and used again.

8.4 Air-in-line Alarm

Cause: When the infusion pump is running, if the air-in-line sensor detects that the size of air bubble is larger than that of the

preset limit, the device will trigger an alarm, stop running, make a high priority alarm sound, the screen appear with the message "Air-in-line" and the operation status indicator flash red at the same time.

Solution:

- 1) Click the " √ " button on the screen to clear the alarm and the message "Air-in-line" disappears.
- 2) To remove air bubbles from the infusion line, close the roller clamp, open pump door, take IV infusion set out, check whether there is air bubbles in the line, shake and move the air bubbles to the drip chamber by hands if there is, reinstall the IV infusion set, close the pump door, open the roller clamp, press the Start/Stop key or click the Start button to restart infusion.
- 3) Check if the air-in-line sensor is clean. If sensor probe is dirty, uninstall the IV infusion set, wipe clean the sensor probe with alcohol, reinstall the IV infusion set, and restart infusion.
- 4) If there is still an Air-in-line alarm, change to a new IV infusion set, install the IV infusion set and restart.
- 5) The air-in-line alarm will be activated too if the infusion line between the infusion bag and the pump is occluded. Remove the occlusion in infusion line and restart infusion.

8.5 Battery Empty Alarm

Cause: When the internal battery is running out, the device will trigger an alarm, stop running, make a high priority alarm sound, the screen appear with the message "Battery Empty" and the operation status indicator flash red at the same time, and the device will stop running and power off after 3 minutes.

Solution: The external power supply should be used immediately. When plugged in the external power supply, the battery charge light goes on and the battery starts charging. When the battery is fully charged, the battery charge indicator goes out.

8.6 Battery & External Power Disconnect Alarm

Cause: When the Infusion pump is running, the external power is disconnected, and the device does not detect the battery signal or the battery is disconnected unexpectedly, the device will trigger an alarm, the screen is black, the operation status indicator flash red at the same time, and the sound and light continue to alarm for 3 minutes before the device automatically power off.

Solution: Use external power supply or battery supply, and restart the device after power supply.

8.7 Motor Error Alarm

Cause: When an error is detected in the motor feedback signal (too slow or too fast, or wrong direction of motor operation etc.), or the sensor detection of push handle position does not match the cumulative amount of infusion, the device will trigger an alarm, stop running, make a high priority alarm sound, the screen appear with the message "Motor Error" and the operation status indicator flash red at the same time.

Solution: Click the "√" button on the screen to clear the alarm. Start the infusion again, still report the fault alarm, please contact our service personnel.

8.8 Communication Error Alarm

Cause: When the communication of the device monitoring CPU is incorrect, the device will trigger an alarm, stop running, make a high priority alarm sound, the screen appear with the message "Communication Error" and the operation status indicator flash red at the same time.

Solution: Click " ✓ " button on the screen to clear the alarm sound. Press the Power on/off key to shut down the device and restart the device. If the fault alarm is still reported, please contact our service personnel.

8.9 Battery Error Alarm

Cause: When the external power is inserted on the device, the device does not detect the battery signal or the battery is disconnected unexpectedly, the device will trigger an alarm, stop running, make a high priority alarm sound, the screen appear with the message "Battery Error" and the operation status indicator flash red at the same time.

Solution: Click " ✓ " button on the screen to clear the alarm sound. Press the Power on/off key to shut down the device and restart the device. If the fault alarm is still reported, please contact our service personnel.

8.10 KVO Alarm

Cause: When KVO is enabled and the infusion is complete, the device will automatically convert to the KVO Rate to continue operation. At the same time, the device will trigger high priority alarm and sound alarm, the screen appear with the message "KVO" and the operation status indicator will flash red at the same

time.

Solution: Click the " ✓ " button on the screen to clear the alarm. The message "KVO" disappear. The device can be reset according to operating steps.

8.11 KVO End Alarm

Cause: When the KVO state runs for 30 minutes until KVO task is completed, the device will trigger an alarm, stop running, make a high priority alarm sound, the screen appear with the message "KVO End" and the operation status indicator flash red at the same time.

Solution: Click the " ✓ " button on the screen to clear the alarm. The message "KVO End" disappear. The device can be reset according to operating steps.

8.12 Standby End Alarm

Cause: When the device is in standby and the standby is over, the device will trigger high priority alarm and sound alarm, the screen appear with the message "Standby End" and the operation status indicator flash red at the same time.

Solution: Click the " ✓ " button on the screen to clear the alarm.

8.13 No AC Power Alarm

Cause: When the device is powered on without the network power and use battery power supply, the device will trigger an alarm, make a low priority alarm sound, the screen appear with the message "No AC Power" and the operation status indicator steady

on yellow at the same time.

Solution: Click the " ✓ " button on the screen or connect to an external power supply to clear the alarm.

8.14 Call Back Alarm

Cause: The device is ready to start the infusion, but the device is not started and the device is placed without operation for the set time, the device will trigger an alarm, make a low priority alarm sound, the screen appear with the message "Call Back" and the operation status indicator steady on yellow at the same time.

Solution: Click the " ✓ " button on the screen to clear the alarm.

8.15 Low Battery Alarm

Cause: When the internal battery is low, the device will trigger an alarm, make a low priority alarm sound, the screen appear with the message "Low Battery" and the operation status indicator steady on yellow at the same time. If the Infusion pump is infusing, the device will not stop infusing.

Solution: The external power supply should be used immediately. When plugged in the external power supply, the battery charge indicator lights up, the battery starts charging, and the message "Low Battery" disappears, the battery icon shows the dynamic effect of charging. When the battery is fully charged, the battery charge light goes out.

8.16 Near End of Infusion Alarm

Cause: When the remaining time is less than or equal to the

set near end of infusion time, the device will trigger an alarm, make a low priority alarm sound, the alarm indicator area at the top of the screen appear with the message "Near End Of Infusion" and the operation status indicator steady on yellow at the same time, and the Infusion pump continues to infusion does not stop.

Solution: Click the "Mute" button on the right of the alarm prompt area to suspend the alarm sound. Focus on the amount of fluid remaining in the Infusion pump and the time remaining, and wait for the Infusion pump infusion to complete.

9 System Parameter Setting

9.1 Bed Number

The Bed number can be set on the Home-Setting- Bed number page.

9.2 Brightness

On the Home- Setting- Brightness page, brightness of the display can be adjusted between 1 and 10 level by clicking on the + or - sign. After brightness setting is completed, click the " ✓ " button to return to the previous page.

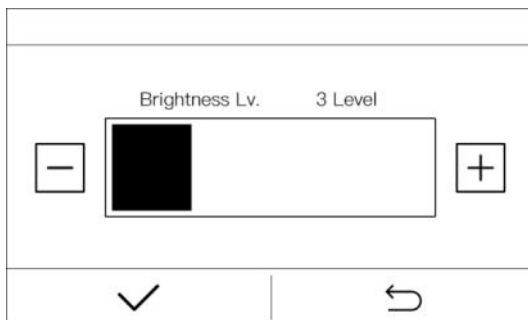


Figure 9-2 Brightness setting

9.3 Alarm Sound Volume

On the Home- Setting- Volume- page, alarm sound volume can be adjusted between 1 and 5 level by clicking on the + or - sign. After sound volume setting is completed, click the " ✓ " button to return to the previous page.

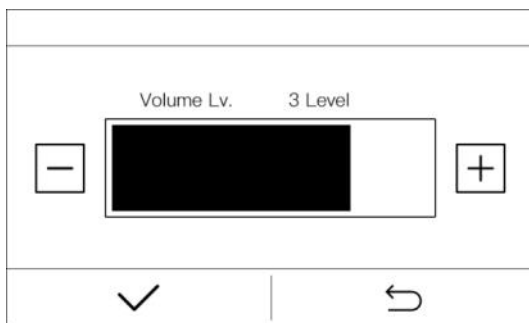


Figure 9-3 Volume setting

9.4 Occlusion Pressure Level

The occlusion pressure Level can be set on the Home- Setting- Level Setting page. Click + or - sign on the screen to choose a pressure level. Click the " ✓ " button to return to the previous page.

Minimum selectable flow rate (0.01mL/h): When the occlusion alarm pressure is at the lowest gear, the alarm time is greater than 5h; When the occlusion alarm pressure is at the highest gear, the alarm time is greater than 20 hours.

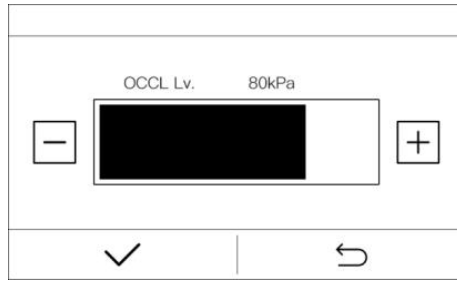


Figure 9-4 Occlusion pressure level setting

9.5 Air Bubble Detection Sensitivity

The Air-in-line detection sensitivity can be set the Home-Setting- Bubble page. Click on the + or - sign to adjust the sensitivity level. After setting is completed, click the " ✓ " button to return to the previous page.

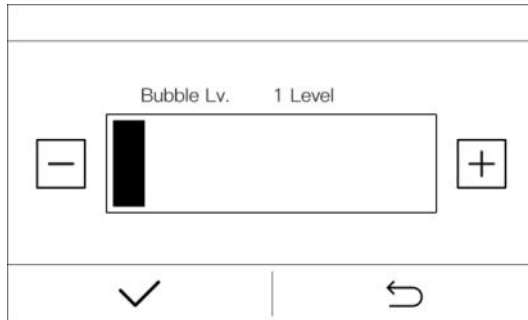


Figure 9-5 Air bubble detection sensitivity setting

The smaller the bubble level, the more sensitive it is. The minimum bubble size detectable for each level is shown in the following table:

Air Bubble Level	1	2	3	4	5	6	7	8	9
Bubble Size (μL)	25	50	100	200	300	500	800	1000	1200

9.6 Bolus Setting

The Bolus Rate and Bolus VTBI can be set on the Home-Setting- Bolus page. Click the corresponding parameter value to set. After the Settings are complete, click the Back button to return to the previous page.

9.7 Purge Setting

The Purge Rate and Purge VTBI can be set on the Home-Setting- Purge page. Click the corresponding parameter value to set. The values set in the Purge setting will not affect the Bolus Rate and Bolus VTBI.

9.8 Call Back Time Setting

On the Home- Setting- Call Back Time page, can set the Call Back Alarm time when the device is placed without operation and not running.

9.9 KVO Setting

On the Home- Setting - KVO page, set the KVO Mode and KVO Rate after the End Of Infusion.

Select the Constant KVO mode and the device will operate according to the currently set KVO Rate.

Select the Variable Speed KVO Mode, the device will automatically determine the KVO Rate level according to the current Infusion rate of the infusion operation, and run according to the KVO Rate of the currently set Infusion rate level.

Constant KVO Mode (V_{KVO} rate setting range: 0.10-5.00 mL/h,

minimum step value 0.01 mL/h)

-When the infusion rate is $>V_{KVO}$: After the infusion task is completed, the pump automatically runs at the V_{KVO} rate to maintain venous patency.

-When the infusion rate is $\leq V_{KVO}$: After the infusion task is completed, the pump continues to run at the infusion rate to maintain venous patency.

Variable Speed KVO Mode (V_{KVO} rate setting range: 0.10-5.00 mL/h, minimum step value 0.01 mL/h)

Users need to set the V_{KVO} rate separately for three scenarios:

- (a) $V_{KVO}(>10)$ rate for infusion rates > 10 mL/h,
- (b) $V_{KVO}(1-10)$ rate for infusion rates between 1-10 mL/h, and
- (c) $V_{KVO}(\leq 1)$ rate for infusion rates ≤ 1 mL/h.

The system will automatically matches the KVO rate based on the infusion rate:

-When the infusion rate is > 10 mL/h: After the infusion task is completed, the pump automatically runs at the $V_{KVO}(>10)$ rate to maintain venous patency.

-When the infusion rate is between 1-10 mL/h: After the infusion task is completed, the pump automatically runs at the $V_{KVO}(1-10)$ rate to maintain venous patency.

-When the infusion rate is ≤ 1 mL/h and the infusion rate is $> V_{KVO}(\leq 1)$ rate: After the infusion task is completed, the pump continues to run at the $V_{KVO}(\leq 1)$ rate to maintain venous patency.

-When the infusion rate is ≤ 1 mL/h and the infusion rate is $\leq V_{KVO}(\leq 1)$ rate: After the infusion task is completed, the pump continues to run at the infusion rate to maintain venous patency.

9.10 Screen Lock Time

On the Home- Setting- Auto Screen Lock page, set the Lock time for the screen and keys of the device during infusion operation.

If the screen lock time is set to 0, the Auto Screen Lock function is disabled.

You can also manually lock the device by clicking the "Lock screen" button on the screen when the device infusion is running. When the screen is locked, click the "Unlock" button on the screen to unlock it. When the screen is locked, other buttons and screen areas are unavailable except the power button and the Unlock button on the screen.

9.11 Near End of Infusion Time Setting

The near end of infusion time can be set on the Home- Setting- NEOI page.

9.12 Prime Prompt Switch

On the Home- Setting- Prime prompt switch page, can choose whether to enable the Prime prompt, if this function is enabled, the "Prime" prompt page will appear when a new infusion task starts after each IV infusion set change.

9.13 History Mode Switch

On the Home- Setting- History mode switch page, can choose whether to enable the History mode, if this function is enabled, the "Sure to load last treatment?" prompt page will be displayed after

each power on.

9.14 IV Set Brand

When the device is not running, click the IV infusion set brand on the upper right corner of the interface, select the corresponding IV infusion set name, and the device will return to the infusion page. The selected brand for the IV infusion set will be shown on the upper right corner of the infusion page, which can remind the operator to use the right IV infusion set to maintain high infusion accuracy.

The calibration operation of IV infusion set, see Section 10.2 Accuracy calibration for IV infusion set.

9.15 Night Mode Setting

On the Setting- Maintenance- Night Mode page, can set the brightness or sound volume for daytime or nighttime, the setting method is the same as that for brightness and sound level.

Night Mode	1/2	
Night Mode Enable		<input checked="" type="checkbox"/>
Volume	3	>
Bright.	3	>
Start Time	19:00:00	>
< > ↶		



Figure 9-15 Night mode setting

9.16 Date/Time Setting

The device of Time and date can be set on the Setting - Maintenance- Date/Time page.

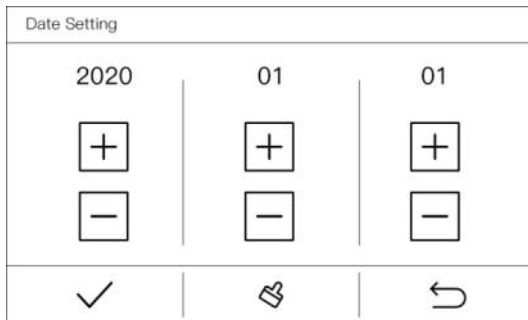


Figure 9-16-1 Date Setting

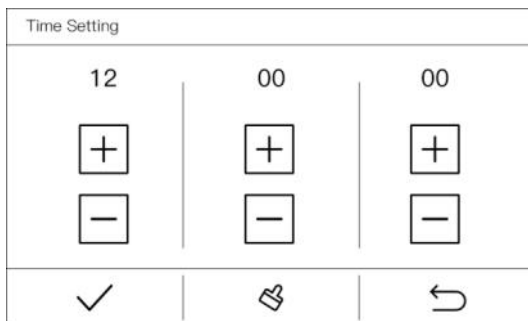


Figure 9-16-2 Time setting

9.17 Maintenance

Maintain the device in "Setting - Maintenance", including the calibration of IV infusion set accuracy, system time settings, system language settings, factory resets, etc. Entering the Maintenance interface requires entering the engineer's maintenance password, password please consult the company's sales and service personnel.

IV infusion set calibration operation instructions see Section 10 of this manual.

10 Accuracy Calibration for IV infusion set

The brand and specification IV infusion sets that have not been calibrated, or IV infusion sets belonging to different production lots with the calibrated IV infusion sets, need to be calibrated before use. In normal use, the IV infusion set should be re-calibrated every 6 months.

10.1 Enter the Infusion set Accuracy Calibration Interface

Enter the system maintenance page by following instructions in "9 System Parameter Setting". Select the brand and model of the infusion set to be calibrated from the Consumables Maintenance-Edit Consumables brand-consumable brand.

10.2 Accuracy Calibration for IV infusion set

- 1) Same as the normal infusion operation, install the IV infusion set first, put the scalp needle into the beaker, place the beaker on the

balance and clear the indicator to zero, and enter the page of "Setting- System Maintenance- Consumable Maintenance - Calibration Consumable".

- 2) Click Brand, select the brand of the IV infusion set to be calibrated, and return to the Calibration Consumable page.
- 3) Click Calibration range, select the high Rate interval or the low Rate (the calibration rate of the low rate is 200mL/h, and the output volume is 10mL; the calibration rate of the high rate is 1000mL/h, and the output volume is 50mL), return to the Calibration Consumable page, and pay attention to the selection of the scalp needle specifications that match the rate.
- 4) Click Calibration Start/Stop, the infusion pump will output a certain amount of solution according to the current calibration interval. When the infusion was completed, the infusion pump automatically stopped running.
- 5) Check the balance reading, convert it to the actual solution volume, and input the actual solution volume in the Volume Output;
- 6) Click the page turning button to the next page, click confirm, and the selection interface will appear. Select " √ " to enter the data saving waiting interface, and wait for the saving to be completed to automatically return to the consumables maintenance page; Select "X" to return to the editing consumables page without saving the current consumables calibration data;
- 7) Return to the next page and click OK to save the current calibration value.

11 Precautions for Using Disposable Consumables

It's suggested to use the recommended consumables. The ambient temperature should be kept at least at 5 °C or above when a recommended consumable is used. The infusion accuracy will be compromised if ambient temperature is lower than 5 °C.

The recommended consumables are listed in the table below:

No.	Brand	Model	Infusion Accuracy	Ambient temperature
1	Hanaco	H-06APD-8	±4.5%	+5 °C ~ +40 °C
2	Wego	Automatic Vent Type With Needle	±4.5%	+5 °C ~ +40 °C
3	Yusheng	ordinary type with needle	±4.5%	+5 °C ~ +40 °C
4	Shinva	ordinary type with needle	±4.5%	+5 °C ~ +40 °C
5	MDK	EF-BS1-P1	±4.5%	+5 °C ~ +40 °C
6	Terumo	1TB*U800B	±4.5%	+5 °C ~ +40 °C

The consumable used must have a medical device product registration certificate, and the consumable specifications are selected in the same specification as the recommended consumable brand. consumable Installation Methods See 7.3.1 Install/replace consumable.

In order to ensure infusion accuracy, when the ambient temperature changes significantly, the device needs to be re-calibrated, calibration method See Section 10 Accuracy Calibration for consumable.


Please strictly follow the requirements described in Section 10 to calibrate and use the consumable when change to a new consumable from a different manufacturer. Otherwise, the infusion accuracy may be compromised.

12 Technical Specification

According to IEC60601-2-24 standard, we tested and obtained the trumpet and flow rate graphs under the conditions of 1mL, 25mL, positive and negative back pressure, etc. The specific data are reflected in the IEC60601-2-24 test reports, and the relevant graphs will update in this section accordingly.


This product is numbered in accordance with the technical requirements of medical device products:

- 1) The methods of controlling Bolus volume before occlusion: The pressure in the occlusion pipeline is released to control the bolus volume by controlling the inversion of the stepper motor.
- 2) Storage time for the electronic memory after power off: same as the product lifetime.
- 3) The maximum volume that the pump can deliver under a single fault condition: 0.4 mL
- 4) Device calibration is measured in ml.
- 5) Minimum flow rate performance curve (1 mL/h)
 - a. The rising curve for HANACO IV infusion set with the minimum flow rate during the first 2 h of the test period.

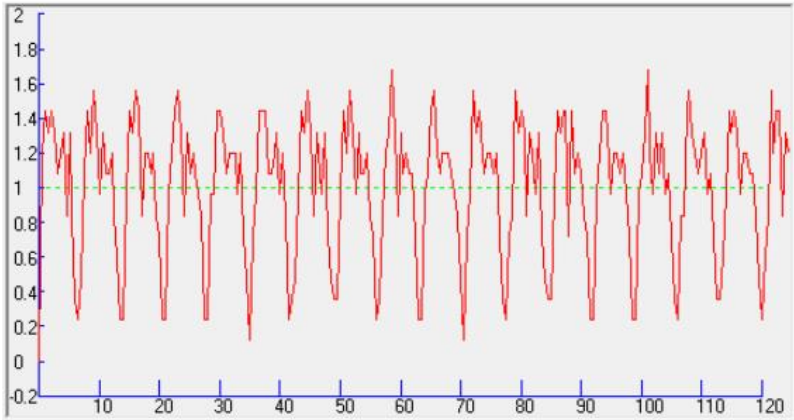

In the above figure, the dashed line shows the set flow rate (1 mL/h in this figure), and the solid line is the continuous connection line for the average flow rate during a sampling period.


- b. The trumpet curve for HANACO IV infusion set with the minimum flow rate during the second hour of the test period.

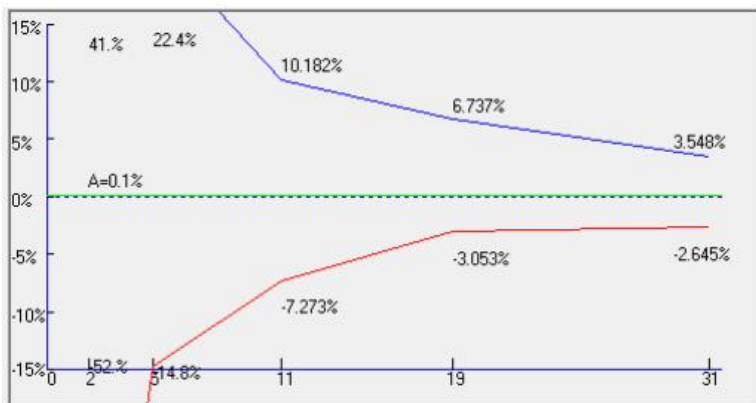

The dashed line in green color is the final value that the infusion error of the device is eventually converging to. The solid blue line above the dashed line is the maximum positive deviation


during the second hour of the test period. The solid red line below the dashed line is the maximum negative deviation during the second hour of the test period.

- c. The trumpet curve for HANACO IV infusion set with the minimum flow rate during the last hour of the run.

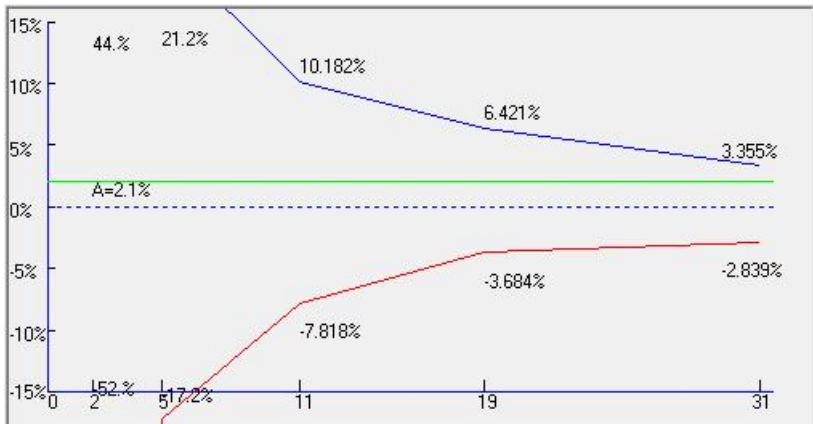


The green dashed line is the value where the device perfusion error finally converges. The solid blue line above the dashed line is the maximum positive deviation within the hour of the run. The solid red line below the dashed line is the maximum negative deviation within the last hour.

- 6) Intermediate flow rate performance curve (25 mL/h)
 - a. The rising curve for HANACO IV infusion set with the intermediate flow rate during the first 2h of the test period.

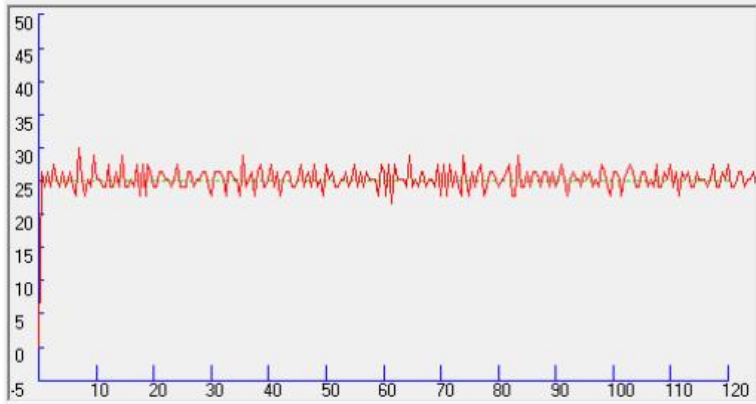

In the above figure, the dashed line shows the set flow rate (25 mL/h in this figure), and the solid line is the continuous connection line for the average flow rate during a sampling period.


- b. The trumpet curve for HANACO IV infusion set with the intermediate flow rate during the second hour of the test period.

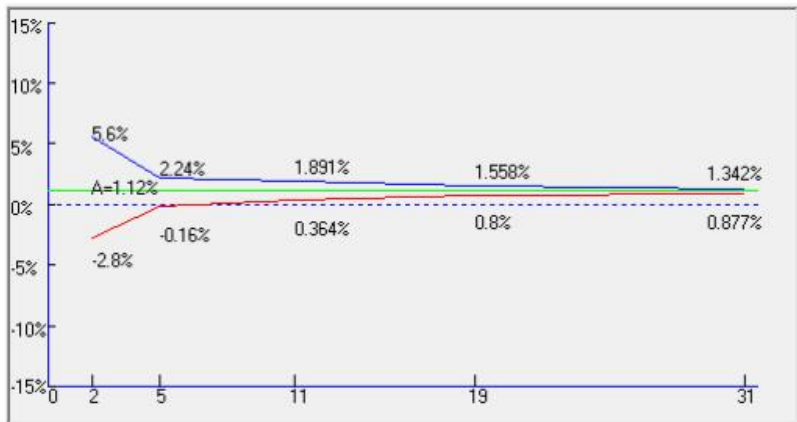

The dashed line in green color is the final value that the infusion error of the device is eventually converging to. The solid blue line above the dashed line is the maximum positive deviation


during the second hour of the test period. The solid red line below the dashed line is the maximum negative deviation during the second hour of the test period.

- c. The trumpet curve for HANACO IV infusion set with the intermediate flow rate during the last hour of the run.

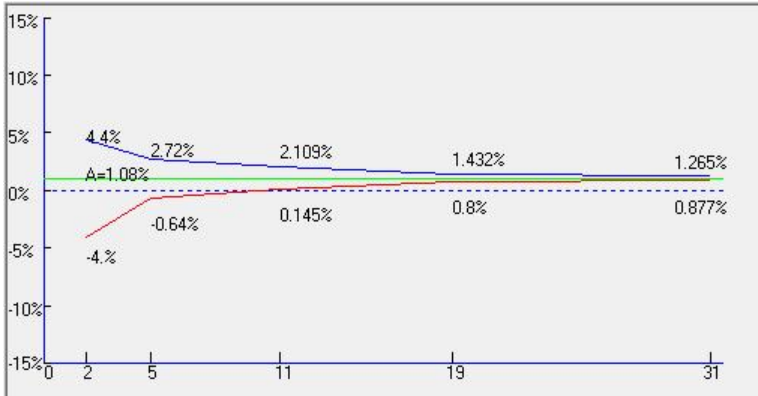


The green dashed line is the value where the device perfusion error finally converges. The solid blue line above the dashed line is the maximum positive deviation within the last hour of the run. The solid red line below the dashed line is the maximum negative deviation within the last hour.

- 7) Intermediate flow rate and back pressure +13.3kPa performance curve
 - a. The rising curve for HANACO IV infusion set with the Intermediate flow rate and back pressure +13.3kPa during the first 2h of the test period.

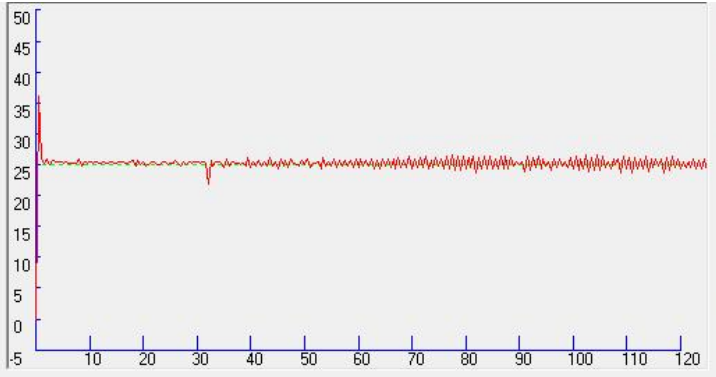

In the above figure, the dashed line shows the set flow rate (25 mL/h in this figure), and the solid line is the continuous connection line for the average flow rate during a sampling period.


- b. The trumpet curve for HANACO IV infusion set with the Intermediate flow rate and back pressure +13.3kPa during the second hour of the test period.

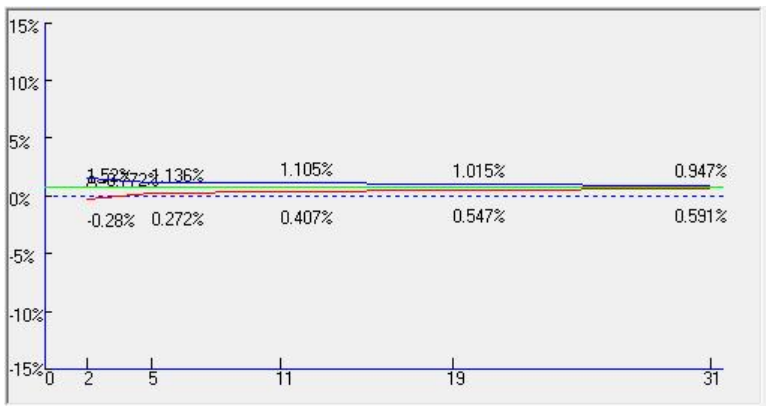

The dashed line in green color is the final value that the infusion error of the device is eventually converging to. The solid blue line above the dashed line is the maximum positive deviation during the second hour of the test period. The solid red line below the dashed line is the


maximum negative deviation during the second hour of the test period.

- 8) Intermediate flow rate and back pressure -13.3kPa performance curve
 - a. The rising curve for HANACO IV infusion set with the Intermediate flow rate and back pressure -13.3kPa during the first 2h of the test period.

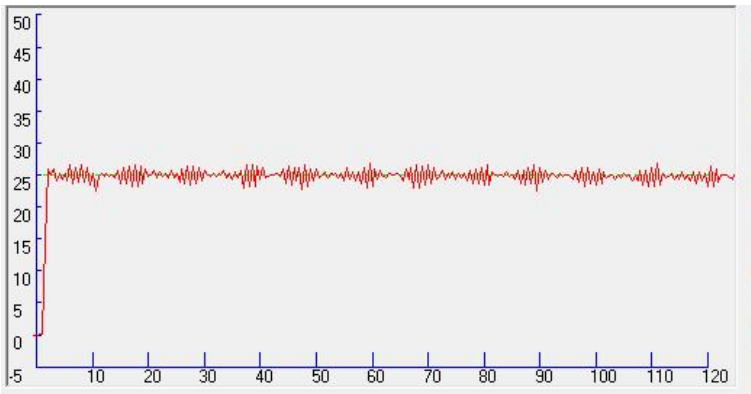


In the above figure, the dashed line shows the set flow rate (25 mL/h in this figure), and the solid line is the continuous connection line for the average flow rate during a sampling period.

- b. The trumpet curve for HANACO IV infusion set with the Intermediate flow rate and back pressure -13.3kPa the second hour of the test period.

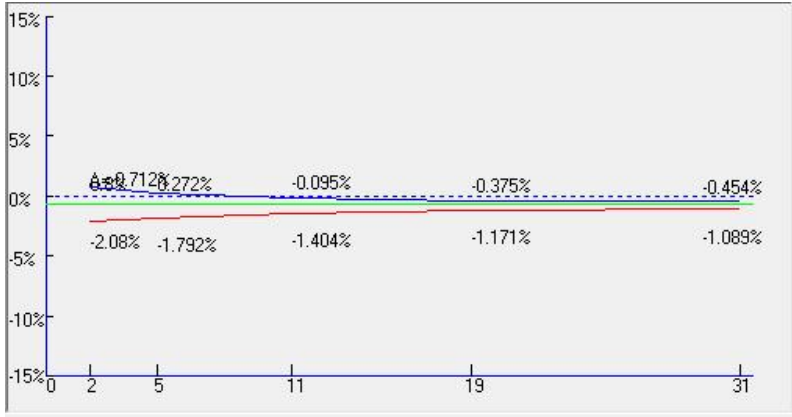


The dashed line in green color is the final value that the infusion error of the device is eventually converging to. The solid blue line above the dashed line is the maximum positive deviation during the second hour of the test period. The solid red line below the dashed line is the maximum negative deviation during the second hour of the test period.

- 9) Intermediate flow rate and below distance 0.5 meter performance curve
 - a. The rising curve for HANACO IV infusion set with the Intermediate flow rate and the supply container below the pump mechanism at a distance of 0.5m during the first 2h of the test period.

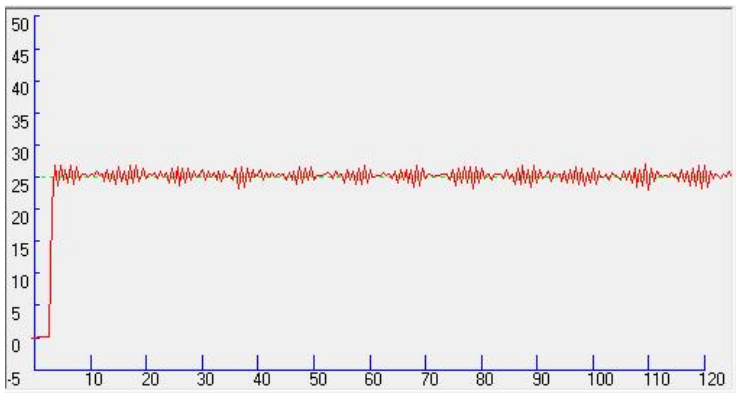

In the above figure, the dashed line shows the set flow rate (25 mL/h in this figure), and the solid line is the continuous connection line for the average flow rate during a sampling period.

- b. The trumpet curve for HANACO IV infusion set with the Intermediate flow rate and the supply container below the pump mechanism at a distance of 0.5m during the second hour of the test period.

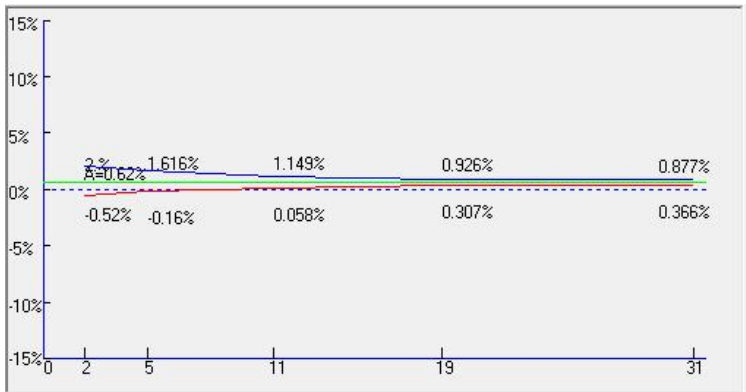

The dashed line in green color is the final value that the infusion error of the device is eventually converging to. The solid blue line above the dashed line is the maximum positive deviation during the second hour of the test period. The solid red line below the dashed line is the maximum negative deviation during the second hour of the test period.

- 10) The sensitivity of the air-in-line sensor: the minimum air bubble size that can be detected is 25 uL.

13 Restore to Factory Setting

Default factory setting as below:

No.	Parameter	Factory presets
1	Brightness level	Level 3

2	System sound level	Level 3
3	Bubble level	Level 3
4	Night mode sound level	Level 3
5	Night mode brightness level	Level 3
6	Occlusion pressure level	120 kPa
7	Night mode	Close
8	Infusion mode	Rate mode
9	KVO	Constant KVO
10	KVO rate	1.00 mL
11	Call Back Time	2 min
12	Near End Of Infusion time	5 min
13	Auto screen lock time	5 min
14	Night mode start time	19:00:00
15	Night mode end time	09:00:00
16	Bolus rate	1200.00 mL/h
17	Purge rate	1200.00 mL/h
18	Bolus volume	5.00 mL
19	Purge volume	15.00 mL

14 Use, Maintenance and Removal of the Internal Battery

The device has an internal rechargeable lithium battery with the following specification: Li-ion Battery ICR18650, capacity: 2200mAh, standard voltage 11.1 V. After the network power is connected, the internal charge management module automatically controls the charging of the lithium battery. When the network power supply is interrupted, the system will automatically switch to the internal battery working state.

When the battery is fully charged and powered by an internal battery, it can operate continuously for more than 6 hours under intermediate rate operating conditions.

Daily maintenance of the battery:

- 1) When the pump is not used for a long time, the internal battery should be fully charged at least once for every 3 months by connecting the device to the mains power to help saving the battery life.
- 2) The device operates at a maximum selectable flow rate (2000.00 mL/h) and can run on batteries for up to 6h10min, making it easy for patient transportation and unexpected power outages.
- 3) Contact the customer service immediately if the internal battery cannot be charged or cannot work normally. Do not disassemble it by yourself. For the healthcare providers who have the ability to repair a device, we will provide training to the related personnel from these facilities.

The device has a internal disposable button battery designed to last longer than 8 years, when the set time is exceeded, need to be disposed with the device in accordance with the instructions for waste disposal in this manual 16.

15 Service and Maintenance

Check the pump before use:

- 1) Check whether there are foreign objects inside the power outlet (such as drug solution residue), and confirm that the device startup self-test is normal.
- 2) Select the correct IV infusion set specification, check the battery power, and charge it in time when the power is low.

During use:

- 1) To avoid giving an incorrect dosage of drug to a patient, please disconnect the pump from the patient before changing a device.
- 2) Please make sure that the infusion line is not kinked. Insert the needle to the vein on a part of the patient' s body where it is not likely to be squeezed or pressed.
- 3) To prevent the spilled drug solution on the pump surface from getting into the inside of the device, wipe it dry immediately if there is a spill.

Storage and daily maintenance:

- 1) To keep the device clean, wipe it clean for at least once a month, which can prevent the corrosion caused by the drug solution and avoid the mobility of the mechanical parts being affected by the dried solution.
- 2) Use a clean and damp cloth or an alcohol pad to wipe clean the device. Take caution to avoid any liquid from entering the device. If disinfection is required, commonly used disinfectants can be used. After using the disinfectant, after wetting with a soft cloth in water, wring out the soft cloth for scrub treatment. When using disinfectants, follow their instructions.
- 3) Check the Low battery Alarm time of the device at least once a month. Make the device standby when the battery is low in non-clinical use, start timing when you hear the alarm of "low battery", and the alarm time should be more than 30 minutes.
- 4) Keep the surface of the air bubble sensor probe clean. A dirty probe will reduce the sensor' s sensitivity in air bubble detection or cause false alarm. Take caution when cleaning the probe to avoid causing any damage.

16 Waste Disposal

16.1 Battery

Please follow local regulations to dispose of used batteries.

16.2 IV infusion set

After use, please dispose of the IV infusion set in accordance with the relevant medical waste disposal regulations.

16.3 Infusion Pump

This device is designed to last 8 years and should be scrapped after it has exceeded the lifetime. End-of-life Infusion pumps can be sent back to the dealer who sold the product or to the Company for proper recycling.

17 Electromagnetic Compatibility

Special precautions regarding Electromagnetic Compatibility (EMC) are required for this device. Must install and use in accordance with the electromagnetic compatibility information specified in this instruction.

Portable and mobile RF communication devices may have an impact on this device.

Must use the cables and accessories provided by this device, and the cable information as follows:

The name of the cable	Length
Power cable	2.9 m

In addition to cables (transducers) sold as spare parts for internal components, the use of accessories and cables (transducers) other than specified may result in an increase in device or system emission or a decrease in immunity.

Devices or systems should not be used close to or stacked with other devices, and if they must be accessed or stacked, observe to verify that they can run normally in the configuration they are using.

The basic performance is to operate on a network power supply (including an internal battery) connection.

Item	Description
The network power supply (including the internal battery) is connected to run	At the intermediate rate of 25.00mL/h and the VTBI of $\geq 10\text{mL}$, start to operation, infusion accuracy error less than $\pm 4.5\%$ and the operation is normal during the process, there should be no abnormal phenomena and failures.

Guidance and manufacture' s declaration - electromagnetic emission

Guidance and manufacture' s declaration - electromagnetic emission

The Infusion pump is intended for use in the electromagnetic environment specified below. The customer or the user of the Infusion pump should assure that it is used in such and environment.

Emission test	Compliance	Electromagnetic environment guidance
RF emissions CISPR 11	Group 1	The device only uses RF energy for its internal functions, so its RF emission is very low and the possibility of interference to nearby electronic devices is very low
RF emissions CISPR 11	Class A	The device is suitable for use in non household and all facilities that are not directly connected to the public low-voltage power supply network of household residential areas.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations /	Not applicable	

flicker emissions IEC 61000-3-3		
------------------------------------	--	--

Guidance and manufacture' s declaration - electromagnetic immunity

Guidance and manufacture' s declaration – electromagnetic immunity

The Infusion pump is intended for use in the electromagnetic environment specified below. The customer or the user of Infusion pump should assure that it is used in such an environment.


Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8KV contact ±15 KV air	±8KV contact ±15KV air	The ground should be made of wood, concrete, or ceramic tiles. If the ground is covered with synthetic materials, the relative humidity should be at least 30%
Electrical fast transient/burst IEC 61000-4-4	±2KV for power supply lines ±1 KV signal input/output	±2 KV for power supply lines Not applicable	The network power supply should have the quality used in typical commercial or hospital environments
Surge IEC 61000-4-5	±1 KV line to line ±2 KV line to ground	±1 KV line to line Not applicable	The network power supply should have the quality used in typical commercial or hospital environments
Voltage dips, short interruptions and voltage variations on power supply	<5%U _T , 0.5 cycle (>U _T , >95% voltage dips) 40%U _T , 5 cycle (> U _T , 60% voltage dips)	<5%U _T , 0.5 cycle (>U _T , >95% voltage dips) 40%U _T , 5 cycle > U _T , 60% voltage dips)	The network power supply should have the quality suitable for typical commercial or hospital environments. If the user of the device needs to operate continuously during

input lines IEC 61000-4-11	70%U _T , 25 cycle (>U _T , 30% voltage dips) <5%U _T 5 s (>U _T , >95% voltage dips)	70%U _T , 25 cycle (>U _T , 30% voltage dips) <5%U _T , 5s (>U _T , >95% voltage dips)	a power outage, it is recommended to use an uninterruptible power supply or battery power supply for the device.
Power frequency (50/60Hz) magnetic field IEC 61000-4-8	400A/m	400 A/m	The power frequency magnetic field should have the horizontal characteristics of the power frequency magnetic field in typical commercial or hospital environments

Note: U_T is the a.c. mains voltage prior to application of the test level.

Guidance and manufacture' s declaration - electromagnetic immunity

Guidance and manufacture' s declaration – electromagnetic immunity			
The Infusion pump is intended for use in the electromagnetic environment specified below. The customer or the user of Infusion pump should assure that it is used in such an environment:			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment guidance
Conduced RF IEC61000-4-6	3 V (Valid value) 150 kHz~ 80MHz (expect the ISM bands between ^a)	10V	Portable and mobile RF communication devices should not be used closer to any part of the device, including cables, than the recommended isolation distance. The distance should be calculated using a formula corresponding to the transmitter frequency.
Radiated RF IEC61000-4-3	10V (Valid value)	10V/m	

	<p>150 kHz~ 80MHz (ISM bands between ^{a)})</p> <p>10V /m 80MHz~2.5GHz</p>		<p>Recommended isolation distance</p> <p>$d = 1.2\sqrt{(P)}$ 150kHz-80MHz $d = 1.2\sqrt{(P)}$ 80MHz-800MHz $d = 2.3\sqrt{(P)}$ 800MHz-2.5GHz</p> <p>Note:</p> <p>P—According to the maximum rated output power of the transmitter provided by the transmitter manufacturer, in watts (W); d —Is the recommended isolation distance in meters (m).^b</p> <p>The field strength of a fixed RF transmitter is determined by surveying the electromagnetic field c, and should be lower than the corresponding level in each frequency range d.</p> <p>Interference may occur near devices marked with the following symbols</p> 
--	--	--	---

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a The ISM bands between 150kHz and 80MHz means 6.765MHz~6.795 MHz、13.553MHz~13.567 MHz、26.957MHz~27.283 MHz 和 40.66MHz~40.70MHz.

b The ISM bands between 150kHz~80MHz and 80M Hz~2.5GHz compliance level, is used to reduce the possibility of interference caused by mobile/portable communication devices accidentally brought into the patient's area. For this purpose, an additional factor of 10/3 is used to calculate the recommended isolation distance of the transmitter within these frequency ranges.

c Fixed transmitter. For example, the field strength of base stations for wireless (cellular/wireless)

telephones and ground mobile radios, business radios, AM and FM radio broadcasts, and television broadcasts cannot be accurately predicted in theory. To evaluate the electromagnetic environment of fixed RF transmitters, consideration should be given to the investigation of electromagnetic sites. If the measured field strength of the location where the infusion pump is located is higher than the applicable RF compliance level mentioned above, the infusion pump should be observed to verify its normal operation. If abnormal performance is observed, additional measures may be necessary, such as adjusting the direction or position of the infusion pump.

d In the entire frequency range of 150kHz to 80MHz, the field strength should be less than 3V/m

Recommended isolation distance between portable and mobile RF communication device and infusion pumps

Recommended isolation distance between portable and mobile RF communication device and infusion pumps

The device is expected to be used in an electromagnetic environment with controlled radio frequency radiation disturbance. Based on the maximum rated output power of communication devices, buyers or users can prevent electromagnetic interference by maintaining the minimum distance between portable and mobile RF communication devices (transmitters) and devices as recommended below.

Maximum rated output power of the transmitter W	Isolation distance corresponding to different frequencies of the transmitter/m			
	150kHz ~ 80MHz expect the ISM bands) $d = 1.2 \sqrt{P}$	150kHz ~ 80MHz (ISM bands) $d = 1.2 \sqrt{P}$	80MHz ~ 800MHz $d = 1.2 \sqrt{P}$	800MHz~2.5GHz $d = 2.3 \sqrt{P}$
0.01	0.12	0.12	0.12	0.23
0.1	0.38	0.38	0.38	0.73
1	1.2	1.2	1.2	2.3
10	3.8	3.8	3.8	7.3
100	12	12	12	23

For the maximum rated output power of transmitters not listed in the above table, it is recommended to isolate the distance d in meters (m), which can be determined by the

formula in the corresponding transmitter frequency column. Here, P is the maximum rated output frequency of the transmitter provided by the transmitter manufacturer, in watts (W).

Note 1: At the frequency points of 80 MHz and 800 MHz, the formula for the higher frequency band is used.

Note 2: The power frequency medical frequency band between 150 kHz and 80 MHz refers to 6.765 MHz to 6795 MHz, 13.553 MHz to 13567 MHz, 26.957 MHz to 27.283 MHz, and 40.66 MHz to 40.70 MHz.

Note 3: The additional factor of 10/3 is used to calculate the recommended isolation distance for transmitters in the engineering and medical frequency bands of 150 kHz to 80 MHz and the frequency range of 80 MHz to 2, 5 GHz, in order to reduce the possibility of interference caused by portable/mobile communication devices accidentally brought into the patient's area.

Note 4: These guidelines may not be suitable for all situations. Electromagnetic propagation is influenced by the absorption and reflection of buildings, objects, and the human body.

18 Antistatic Precautions

The Infusion pumps have been tested and comply with medical device standard IEC 60601-1-8.

When using this device, the user should not touch the pins of connectors marked with an electrostatic discharge warning symbol and should not connect to these connectors unless electrostatic discharge precautions are used.

The operator should be aware of the following things:

a) Unless appropriate preventive measures have already been taken, do not use hand or hand tool to touch connectors with electrostatic discharge warning signs. Preventive measures include: 1) Methods for preventing electrostatic charge accumulation (such as air conditioning, air humidification, floor conductive coating or Non-synthetic clothing); 2) Discharge electrostatic charge from human body to the framework of device, or to the ground, or to a

large piece of metal; 3) Use a wrist band to connect human body to the device or to the ground.

b) All staff who may be in contact with connectors with electrostatic discharge warning signs should receive training, including all clinical/biomedical engineering and healthcare personnel.

c) Electrostatic discharge training should include the introduction of static charges in the theory of physics, the voltage that may be produced in normal practice, and the damage to the electronic components caused by the electrostatic charge from an operator. Further, methods for how to prevent electrostatic charge accumulation should be provided, as well as how and why to discharge the electrostatic from human body to the framework of device or to the ground, and how to use wrist band to connect someone' s body to the device or to the ground.

19 Packaging and Accessories

The list of recommended accessories for use with this device (single unit) is as follows:

Accessories	Qty	Unit
IFU	1	book
Power cable	1	set
Fixation clamp	1	pcs
Other accessories can be found in the packing list.		



Legal manufacturer: MDKMed Medical Technology Co., Ltd.

Address: 502A, Building 7, No. 22, Xinyan Road, Donghu Street,
Linping District, Hangzhou City, Zhejiang Province, P. R China 311323

After sale service: MDKMed Medical Technology Co., Ltd.

Tel: 400-880-8392

Email: sales@graseby.com

Web: <http://www.graseby.com>

Revision date: 02-01-2025

Part No.: 07.01.0168A



European Authorized Representative:

MedNet EC-REP C IIb GmbH

Borkstrasse 10

48163 Münster, Germany